

Code No. 27402

## Human Total HGF Assay Kit - MCM

### INTRODUCTION

Hepatocyte growth factor/scatter factor (HGF) is a liver regeneration and growth factor that was isolated and purified from the plasma of the person who suffers fulminant hepatitis. In addition blood and tissue concentrations being assessed during liver injury, increases in concentration in body fluids have been reported in other types of disease besides liver disease (e.g., inflammatory diseases, neoplastic diseases, and fibrosis). HGF is secreted as an inactive single-chain pro-form, and it exerts its physiological activity after processing by HGF activator (HGFA), coagulation factor VIIa, etc., and conversion to activated HGF, a heterodimer. This product is a kit that measures human HGF by the WHO standard.

### PRINCIPLE

This kit is a solid phase sandwich ELISA using 2 kinds of high specific antibodies. Tetra Methyl Benzidine (TMB) is used as coloring agent (Chromogen). The strength of coloring is in proportion to the quantities of human total HGF.

### MEASUREMENT RANGE

78.13 ~ 5,000 pg/mL

### INTENDED USE

For research use only, not for use in diagnostic procedures.

- The Human Total HGF Assay Kit is a complete kit for the quantitative determination of human total HGF in serum, EDTA-plasma or cell culture media.
- Determination of human total HGF is affected by the presence of heparin in samples, so please use EDTA-plasma as a sample instead of heparin plasma.

### KIT COMPONENT

1	Precoated plate	: Anti- Human Total HGF Rabbit IgG Affinity Purify	96Well x 1
2	Labeled antibody Conc.	: (30X) HRP conjugated Anti-Human HGF Goat IgG Fab' Affinity Purify	0.4mL x 1
3	Standard	: Human Total HGF	0.5mL x 2
4	EIA buffer	: 1% BSA, 0.05% Tween 20 in PBS	30mL x 1
5	Solution for Labeled antibody:	1% BSA, 0.05% Tween 20 in PBS	12mL x 1
6	Chromogen	: TMB solution	15mL x 1
7	Stop solution	: 1N H <sub>2</sub> SO <sub>4</sub>	12mL x 1
8	Wash buffer Conc.	: (40X) 0.05% Tween20 in phosphate buffer	50mL x 1

### OPERATION MANUAL

#### 1. Materials needed but not supplied

- Plate reader (450nm)
- Graduated cylinder and beaker
- Incubator (37°C±1°C)
- Paper towel
- Washing bottle for precoated plate
- Disposable test tube for "2, Labeled antibody Conc." and "6, Chromogen"
- Micropipette and tip
- Deionized water
- Graph paper (log/log)
- Tube for dilution of Standard


#### 2. Preparation

- 1) Preparation of wash buffer  
"8, Wash buffer Conc." is a concentrated (40X) buffer. The temperature of "8, Wash buffer Conc." shall be adjusted to room temperature and then, mix it gently and completely before use. Dilute 50 mL of "8, Wash buffer Conc." with 1,950mL of deionized water and mix it. This is the wash buffer for use. This prepared wash buffer shall be stored in refrigerator and used within 2 weeks after dilution.
- 2) Preparation of Labeled antibody  
"2, Labeled antibody Conc." is a concentrated (30X). Dilute "2, Labeled antibody Conc." with "5, Solution for Labeled antibody" in 30 times according to required quantity into a disposable test tube. Use this resulting solution as Labeled antibody.  
Example)  
In case you use one slit (8 well), the required quantity of Labeled antibody is 800 µL. (Dilute 30 µL of "2, Labeled antibody Conc." with 870 µL of "5, Solution for Labeled antibody" and mix it. And use the resulting solution by 100 µL in each well.)  
This operation should be done just before the application of Labeled antibody.  
The remaining "2, Labeled antibody Conc." should be stored at 4°C in firmly sealed vial.
- 3) Preparation of Standard  
Put just 0.5 mL of deionized water into the vial of "3, Standard" and mix it gently and completely. This solution is 10,000 pg/mL Human Total HGF standard.
- 4) Dilution of Standard  
Prepare 8 tubes for dilution of "3, Standard". Put 230µL each of "4, EIA buffer" into the tube.  
Specify the following concentration of each tube.

Tube-1	5,000 pg/mL
Tube-2	2,500 pg/mL
Tube-3	1,250 pg/mL
Tube-4	625 pg/mL
Tube-5	312.5 pg/mL
Tube-6	156.25 pg/mL
Tube-7	78.13 pg/mL
Tube-8	0 pg/mL (Test Sample Blank)

Put 230 µL of Standard solution into tube-1 and mix it gently. Then, put 230 µL of tube-1 mixture into tube-2. Dilute two times standard solution in series to set up 7 points of diluted standard between 5,000 pg/mL and 78.13 pg/mL. Tube-8 is the test sample blank as 0 pg/mL.

See following picture.

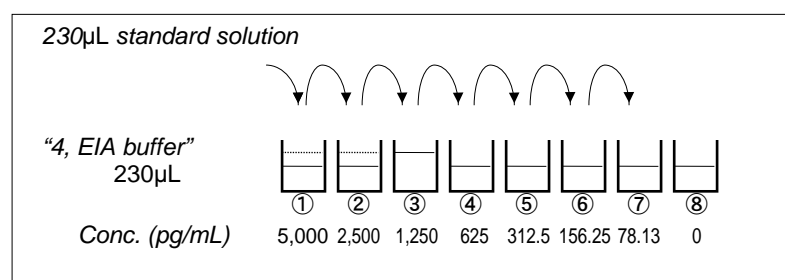

#### 5) Dilution of test sample

Test sample should be diluted with "4, EIA buffer" as the need arises.

If the concentration of Human Total HGF in samples may not be estimated in advance, the pre-assay with several different dilutions will be recommended to determine the proper dilution of samples.

### 3. Measurement procedure

All reagents shall be brought to room temperature approximately 30 minutes before use. Then mix it gently and completely before use. Make sure of no change in quality of the reagents. Standard curve shall be prepared simultaneously with the measurement of test samples.

Reagents	Test Sample	Standard	Test Sample Blank	Reagent Blank
	Test sample 100 µL	Diluted standard (Tube 1~7) 100 µL	EIA buffer (Tube-8) 100 µL	EIA buffer 100 µL
Incubation for 60 minutes at 37°C with plate lid				
Washing 7 times				
Labeled Antibody	100 µL	100 µL	100 µL	-
Incubation for 30 minutes at 37°C with plate lid				
Washing 9 times				
Chromogen	100 µL	100 µL	100 µL	100 µL
Incubation for 30 minutes at room temperature (shielded)				
Stop solution	100 µL	100 µL	100 µL	100 µL
Read the plate at 450nm against a Reagent Blank within 30 minutes after addition of Stop solution.				

- 1) Determine wells for reagent blank. Put 100 µL each of "4, EIA buffer" into the wells.
- 2) Determine wells for test sample blank, test sample and diluted standard. Then, put 100 µL each of test sample blank (tube-8), test sample and dilutions of standard (tube-1~7) into the appropriate wells.
- 3) Incubate the precoated plate for 60 minutes at 37°C after covering it with plate lid.
- 4) Wash each well of the precoated plate vigorously with wash buffer using the washing bottle. Then, fill each well with wash buffer and leave the precoated plate laid for 15~30 seconds. Remove wash buffer completely from the precoated plate by snapping. This procedure must be repeated more than 7 times. Then, remove the remaining liquid from all wells completely by snapping the precoated plate onto paper towel.  
*In case of using a plate washer, after 4 times washing with plate washer, washing with above washing bottle must be repeated 3 times.*
- 5) Pipette 100 µL of labeled antibody solution into the wells of test samples, diluted standard and test sample blank.
- 6) Incubate the precoated plate for 30 minutes at 37°C after covering it with plate lid.
- 7) Wash the precoated plate 9 times in the same manner as 4).
- 8) Taken the required quantity of "6, Chromogen" into a disposable test tube. Then, pipette 100 µL from the test tube into the wells. Please do not return the rest of test tube to "6, Chromogen" bottle to avoid contamination.
- 9) Incubate the precoated plate for 30 minutes at room temperature in the dark. The liquid will turn blue by the addition of "6, Chromogen".
- 10) Pipette 100 µL of "7, Stop solution" into the wells. Mix the liquid by tapping the side of precoated plate. The liquid will turn yellow by the addition of "7, Stop solution".
- 11) Remove any dirt or drop of water on the bottom of the precoated plate and confirm there is no bubble on the surface of the liquid. Then, run the plate reader and conduct measurement at 450nm against a reagent blank. The measurement shall be done within 30minutes after the addition of "7, Stop solution".

### SPECIAL ATTENTION

- 1) Test samples should be measured soon after collection. For the storage of test samples, store them frozen and do not repeat freeze/thaw cycles. Thaw the test samples at low temperature and mix them completely before measurement.
- 2) Test samples should be diluted with "4, EIA buffer", if the need arises.
- 3) Duplicate measurement of test samples and standard is recommended.
- 4) Use test samples in neutral pH range. The contaminations of organic solvent may affect the measurement.
- 5) Use only wash buffer contained in this kit for washing the precoated plate. Insufficient washing may lead to the failure in measurement.


- 6) Remove the wash buffer completely by tapping the precoated plate on paper towel.  
Do not wipe wells with paper towel.
- 7) "6, Chromogen" should be stored in the dark due to its sensitivity against light.  
"6, Chromogen" should be avoided contact with metals.
- 8) Measurement should be done within 30 minutes after addition of "7, Stop solution".

#### CALCULATION OF TEST RESULT

Subtract the absorbance of test sample blank from all data, including standards and unknown samples before plotting. Plot the subtracted absorbance of the standards against the standard concentration on log-log graph paper. Draw the best smooth curve through these points to construct the standard curve. Read the concentration for unknown samples from the standard curve.

Example of standard curve

Conc. (pg/mL)	Absorbance (450nm)
5,000	2.370
2,500	1.216
1,250	0.610
625	0.323
312.5	0.178
156.25	0.115
78.13	0.081
0 (Test Sample Blank)	0.045

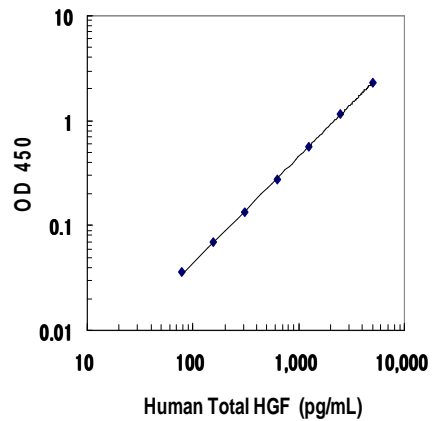

\* The typical standard curve is shown above. This curve can not be used to derive test results. Please run a standard curve for each assay.

#### PERFORMANCE CHARACTERISTICS

##### 1. Titer Assay (Samples with standard added are used.)

Specimen	Titer (X)	Measurement Value (pg/mL)	Theoretical Value (pg/mL)	%
10 % FCS added RPMI-1640	4	2,554.56	2,533.85	100.8
	8	1,198.83	1,254.44	95.6
	16	582.44	641.38	90.8
Human Serum	2	2,203.12	2,668.44	82.6
	4	1,238.30	1,332.81	92.9
	8	690.08	671.90	102.7
Human Plasma (EDTA)	2	2,567.82	2,572.87	99.8
	4	1,215.87	1,287.01	94.5
	8	560.30	649.81	86.2

##### 2. Added Recovery Assay

Specimen	Theoretical Value (pg/mL)	Measurement Value (pg/mL)	%
10 % FCS added RPMI-1640 (x4)	2,500.00	2,441.02	97.6
	1,250.00	1,181.96	94.6
	625.00	595.66	95.3
Human Serum (x2)	449.36	404.08	89.9
	293.11	257.43	87.8
	214.99	201.51	93.7
Human Plasma (EDTA) (x2)	2,557.20	2,461.04	96.2
	1,307.20	1,229.84	94.1
	682.20	564.38	82.7

##### 3. Intra - Assay

Measurement Value (pg/mL)	SD value	CV value (%)	n
1,625.73	111.74	6.9	24
843.27	58.01	6.9	24
226.64	19.87	8.8	24

##### 4. Inter - Assay

Measurement Value (pg/mL)	SD value	CV value (%)	n
1,570.10	157.44	10.0	32
831.69	60.44	7.3	32
221.74	20.44	9.2	32

#### 5. Specificity

Compound	Cross Reactivity
Human Total HGF	100.0%
Human HAI-1	≤0.1%
Human HAI-2	≤0.1%
Human c-Met	≤0.1%
Human HGFA	≤0.1%

#### 6. Sensitivity

11.00 pg/mL

The sensitivity for this kit was determined using the guidelines under the National Committee for Clinical Laboratory Standards (NCCLS) Evaluation Protocols. (National Committee for Clinical Laboratory Standards Evaluation Protocols, SC1, (1989) Villanova, PA: NCCLS.)

#### PRECAUTION FOR INTENDED USE AND/OR HANDLING

1. All reagents should be stored at 2-8°C. All reagents shall be brought to room temperature approximately 30 minutes before use.
2. "3, Standard" is lyophilized products. Be careful to open this vial.
3. "7, Stop solution" is a strong acid substance. Therefore, be careful not to contact your skin and clothes with "7, Stop solution" and pay attention to the disposal of "7, Stop solution".
4. Dispose used materials after rinsing them with large quantity of water.
5. The precipitation may grow in "2, Labeled antibody Conc.", however, there is no problem in the performance.
6. Wash hands after handling reagents.
7. Do not mix the reagents with the reagents from different lot or different kit.
8. Do not use the reagents expired.
9. This kit is for research purpose only. Do not use for clinical diagnosis.

#### STORAGE AND THE TERM OF VALIDITY

Storage Condition : 2 ~ 8°C

The expiry date is specified on outer box.

#### REFERENCE

1. Kataoka H, Hamasuna R, Itoh H, Kitamura N, Kono M. Activation of hepatocyte growth factor/scatter factor in colorectal carcinoma. *Cancer Res.* 2000 Nov 1;60(21):6148-59.
2. Shimomura T, Denda K, Kitamura A, Kawaguchi T, Kito M, Kondo J, Kagaya S, Qin L, Takata H, Miyazawa K, Kitamura N. Hepatocyte growth factor activator inhibitor, a novel Kunitz-type serine protease inhibitor. *J Biol Chem.* 1997 Mar 7;272(10):6370-6.
3. Kawaguchi T, Qin L, Shimomura T, Kondo J, Matsumoto K, Denda K, Kitamura N. Purification and cloning of hepatocyte growth factor activator inhibitor type 2, a Kunitz-type serine protease inhibitor. *J Biol Chem.* 1997 Oct 31;272(44):27558-64.

Version 1.5

Made in Japan.



Distributed By:

**IBL-America, Inc.**

8201 Central Ave NE, Suite P

Minneapolis, MN 55432, USA

[info@ibl-america.com](mailto:info@ibl-america.com)

(888) 523 1246