

Manufactured for Immuno-Biological Laboratories Inc.
(IBL-America)
8201 Central Avenue, NE, Suite P
Minneapolis, MN 55432
Tel: 763-780-2955
Toll Free: 1-888-523-1246



Adiponectin


ELISA

Enzyme Immunoassay for determination of
human Adiponectin
English

For Research Use Only.
Not for use in diagnostic procedures.



Gesellschaft für Forschung und Herstellung von Diagnostika GmbH

 : Aspenhastr. 25 • D-72770 Reutlingen / Germany
Phone: + 49 - (0) 7121 51484-0 • Fax: + 49 - (0) 7121 51484-10
E-mail: contact@mediagnost.de • <http://www.mediagnost.de>

Symbols/ Symbole /Symboles/ Simboli/ Símbolos/ Símbolos/ symbolen/ Symboler/ Symboler/Symbole/ Szimbólumok/ Symboly/ Symboly/ Символы/ Sümbolid/ Σύμβολα/ Simboluri/ Simboli/ Symbolit

DIN EN ISO 15223-1



Expiry date/ Verfallsdatum/ Date de péremption/ Data di scadenza/ Fecha de caducidad/ Data de validade/ Uiterste gebruiksdatum/ Udløbsdato/ Bäst före-datum/ Termin ważności/ Lejárati idő/ Čas expirácie/ Doba expirace/ Срок на годност/ Αεγυμίσκυρπέεβ/ Ημερομηνία λήξης/ Data de expirare/ Rok uporabe/ Viimeinen käyttöpäivä



Consider instructions for use/ Bitte Gebrauchsanweisung beachten/ Consultez la notice d'utilisation/ Consultare le istruzioni per l'uso/ Consulte las instrucciones de uso/ Respeitar as instruções de utilização./ A.u.b. de gebruiksaanwijzing volgen/ Se brugsanvisningen/ Läs anvisningarna före användning/ Proszę przeczytać instrukcję obsługi/ Vegye figyelembe a használati utasításban foglaltakat/ Postupujte podľa pokynov na použitie/ Dodržujte návod k použití/ Моля, спазвайте инструкцията за употреба/ Palun järgige kasutusjuhendit./ Λάβετε υπόψη σας τις οδηγίες χρήσης/ Vá rugám sã respectați instrucțiunile de utilizare/ Upoštečajte navodila za uporabo! Lue käyttöohje huolellisesti!

gg

Lot-Batch Number/ Charge-Chargennummer/ Lot-Code du lot/ Lotto-Numero di lotto/ Lote-Código de lote/ Lote-Código do Lote/Lot-Partijnummer/ Lot-Batchkode/ Partisatskod/ Numer serii/ Tétel-sarzs szám/ Číslo šarže/ Číslo šarže/Партиден номер/Partii – partii number/ Παρτίδα-αριθμός παρτίδας/ Lot-număr lot/ Številka serije/ erä



Manufactured by/ Hergestellt von/Fabriqué par/ Prodotto da/Fabricado por/ Fabricado por/ Vervaardigd doo/Fabrikation af /Tillverkad av/ Wyprodukowane przez/ Gyártotta / Vyrobené/ Vyrobeno v/ Производител/ Tootja/Κατασκευάζεται από/ Produs de/Proizvajalec/ Valmistaja

h

Catalogue Number/ Bestellnummer/ Numéro de référence/Numero di riferimento/Número de referencia/ Número de Referência/ Referentienummer/Referenzenummer /Beställningsnummer/ Numer katalogowy/ Rendelési szám/Katalógové číslo/ Objednací číslo/Καταλογος номер/Tellimnumber/ Αρ. παραγγελίας/Număr comandă/ Številka naročila/ viite tai tilausnumero

l

Store at between/ Lagerung bei zwischen/ Conserver à entre/ Conservare a tra/ Conservar a temp. Entre/ Armazenaer entre/ Bewaar bij tussen/ Opbevares mellem/ Förvaras vid/ Przechowywać w/ Tárolási tartomány/ Skladujte v rozsahu / Skladujte v rozmezi/ Температурно ограничение/ Säilitada temperatuuridel/ Φύλαξη σε θερμοκρασία/ Depozitare între/ Skladiščenje med/ Säilytys x-y Celsiusasteen lämpötilassa

X

Contains sufficient for x tests/ Inhalt ausreichend für x Tests/ Contient suffisant pour x tests/ Contenuto sufficiente per x test/ Contenido suficiente para x pruebas/ Conteúdo suficiente para x testes/ Bevat voldoende voor x bepalingen/ Indeholder tilstrækkeligt til x prøver/Innehållet räcker till x analyser/ Zawartość na x testów/ Tartalma x tesztt elvégzésére elegendő/ Obsahuje materiál pre x testov / Obsah dostacuje pro x testů/ Съдържание достатъчно за x тестове/ Sisust jätkub x katse jaoks/ Το περιεχόμενο επαρκεί για x δοκιμές/ Conținut suficient pentru x teste/Vsebina zadostuje za x preizkusov/ Sisältö riittää x testille



Keep away from sunlight/ Nicht dem Sonnenlicht aussetzen/ Conserver à l'abri de la lumière/ Conservare al riparo della luce solare/ No exponer a la luz solar/ Proteger da luz solar/ Niet aan zonlicht blootstellen/ Må ikke udsættes for sollys/ Utsätt inte för solljus/ Nie wystawiać na słońce/ Napfénytől távol tartandó/ Nevystavovat slnečnému svetlu/ Nevystavovat slunečnímu světlu/ Да се предпазва от слънчева светлина/ Kaista otsese päikesekiirguse eest/ Κρατήστε το μακριά από την ηλιακή ακτινοβολία/ ΤΙνετj departe de lumina soarelui/ Ne izpostavljajte sončni svetlobi/ suojaa auringonvalolta



Incubation time/ Inkubationszeit/ Temps d'incubation/ Tempo d'incubazione/ Tiempo de incubación/ Tempo de incubação/ incubatietijd/Inkubationstid/ inkubationstid/ Czas inkubacji/ Inkubációs idő/ Inkubačná lehota/ Inkubační doba/ Инкубационен период/ Inkubatsiooniaeg/ Χρόνος επώασης/ Timp de incubare/Inkubacijska doba/ inkubaatioaika

°C

incubate at / Inkubation bei/ Incuber à/Incubare a/incubar a/Incubar a/incubatiemtemperatuur/Inkubation ved/inkubacja vid/Inkubacja przy/Inkubáció hőmérséklete/Inkubácia pri/Inkubace při/Инкубира се при/Inkubatsioon temperatuuril/Επώαση στους/Incubare la/Inkubacija pri/ inkubaatiolämpötila



Shaking/ Schütteln/ Mélanger/ Agitare/ Agitar/ Agitação/ Schudden/ Ryster/ Skaka/ Wstrząsanie/ Rázás/ Pretrepat/ Protřepat/ Разклацане/ Raputada/ Ανακινήστε/ Vibrare/ Stresite/ Sekoita

MTP

Mikrotiterplate/ Mikrotiterplatte/ plaque de microtitrage/ Piastra di microtitolazione/ Placa de microtitulación/ Placa de Microtitulação/ microtiterplaat/ Mikrotiterplade/ mikrotiterplatta/ microtiterplaat/ Płytko microtiter/ Mikrotiter lap/ Mikrotitračná podložka/ Mikrotitrační podložka/ Микротитърна плака/ Mikrotiiterplaat/ Τρυβλίο μικροπιλοδότησης/ Microplacã/ Mikrotitrska plošca/ Mikrotittrauslevy

Rec in

Reconstitute in/ Rekonstituieren in/ Reconstituer dans/ Ricostituire ne/ Reconstituir en/ Reconstituir em/ rekonstituieren in/ Rekonstituér i/ rekonstituera/ Rekonstytuować w/ Helyeállitás/ Znovu připravit za/ Znovu připravit za/ Разтваряне в/ Moodustada uesti / Ανασυστήστε σε/ Reconstituire în/ Predelava v /rekonstituoi

SPE

Sample/ Probe /Echantillon/ campione/ Muestra/ Amostra/ monster/ Prové/ prov/ Próbka/ Minta/ Vzorka/ Vzorek/ Проба/ Proov/ Δείγμα/ Probã/Vzorec/Näyte

AbCONJ

AK Antibody and Enzyme Conjugate/ Antikörper und Enzym Konjugat/ anticorps conjuguée et conjugué enzymatique/ Coniugato di anticorpo ed enzima/ Conjugado de anticuerpos y enzimas/ Conjugado Anticorpo-Enzima/ antilichaam- en enzymconjugaat/ Antistoffer og enzym-konjugat/ antikropp- och enzymkonjugat (antikropp och enzym, konjugat)/ Konjugat antycial i enzymów/ Antitest és enzim páros/ Protílátkový a enzymatický konjugát/ Protílátkový a enzymatický konjugát/ Антитяло и ензим конюгат/ Antikehad ja ensüümi konjugaat/ Σύμπλοκο αντισώματος-ενζύμου/ Compuși din anticorpi și enzime/ Antitelesa in konjugat encima/ Vasta-aine ja entsymi konjugaatti

BUF X

VP Buffer / puffer/ Tampon / Tampone / Tampón / Tampão / buffer/ buffer/ buffert/ Bufor / puffer/ pufer/ pufr/ буфер/ puhver/ puhmistικό διάλυμα / Tampon /pufer/ puskuriliuos

DILU BUF X

Dilute in Buffer X/ Verdünnen in Puffer X/ Diluer dans le tampon X/ Diluire nel tampone X/ Diluir en tampón X/ Diluir no Tampão X/ verdunnen in buffer X/ Fortyndes i buffer X/ spädi i buffert X/ Rozcieńczenie w buforze X/ Hígítás X pufferben/ Riedit i pufrü X/ Ředit v pufru X/ Разреждане в буфер X/ Lahjendada puhvrís X/ Αραιώστε σε puhmistικό διάλυμα X/ Diluați în tamponul X/ Razredčiti v pufru X/ Iaimennetaan x puskuriin

STD X

A-E Standard X /Standard X/ Etalon X/ Standard X/ Estándar X/ Standard X/ standaard X/ Standard X/ standard X/ Standard X/ Standard X/ Štandard X/ Standard X/ Стандарт X/ Standard X/ Πρότυπο X/ Standard X/ Standardni X/ Kalibraattori X

Control	KS1/ KS2	Control Serum / Kontrollserum/ Contôle sérique/ Siero di controllo/ Suero de control/ Soro de Controllo/ controleserum/ Kontrolserum/ Kontrollserum/ Serum kontrolne/ Ellenőrző szérum/ Kontrolné sérum/ Kontrolní sérum/ Контролен серум/ Kontrollseerum/ Ορός ελέγχου/Ser de control/ Kontrolni serum/ Kontrolli seerumi
WASHBUF 20x	WP	Washing Buffer Concentrate/ Waschpufferkonzentrat/ Tampon de lavage conc./ Tampone di lavaggio concentrato/ Tampón de lavado concentrado/ Tampão de Lavagem Concentrado/ wasbuffer, geconcentreerd/ Vaskebufferkonzentrat/ Vaskebufferkonzentrat/ tvättbuffertkonzentrat/ Bufor płukania koncentrat/ Mosópuffer koncentrátum/ Концентрат вымывающего пuffers/ Концентрат на промивен бѹфер/ Pesupuhvi kontsentraat/ Συμπύκνωμα ρυθμιστικού διαλύματος έκπλυσης/ Concentrat pentru tamponul de spălare/ Koncentrat izpiralnega puffer/ Pesuliuositiiviste
WASHBUF		Washing Buffer / Waschpuffer/ Tampon de lavage/ Tampone di lavaggio/ Tampón de lavado/Tampão de Lavagem/ wasbuffer/ Vaskebuffer/ tvättbuffert/ Bufor płukania/ Mosópuffer/ Vymývací puffer/ Vymývací pufr/ Промивен бѹфер/ Pesupuhver/ Ρυθμιστικό διάλυμα έκπλυσης/ Tampon pentru spălare /Izpiralni puffer/ Pesuliuos
SUBST TMB	S	Substrate/ Substrat/ Substrat/ Substrato/ Substrato/ Substrato/ substraat/ Substrat/ Substrat/ Substrat/ Szubsztrátum/ Substrát/ Substrát/ Субстрат Substraat/ Υπόστρωμα/ Substrat/ Substrat/ Substraatilius
H₂SO₄	SL	Stop Solution/ Stopp Lösung/ Stop Solution/ Soluzione di stop/ Stop Solución/ Solução Stop/ stopoplossing/ Stopopløsning/ Stopplösning/ Stop roztwór/ Megállító oldat/ Roztok na ukončenie/ Roztok pro ukončení/ Стопирещ разтвор/ Stopp-lahus/ Διάλυμα διακοπής/ Soluție de oprire/ Stop raztopina/ Pysäytysliuos
TAPE		Cover Plate with sealing tape /Platte abkleben/ Recouvrir la microplaque avec bande adhésive/ Coprire la piastra con nastro adesivo/ Cubrir la placa con una cinta adhesiva/ Cobrir a Placa com fita adesiva/ plaatje met tape afdekken/ Afdækningsplade med tape/ maskera platta/ Odkleić płytke/ Tányér leragasztása/ Oblepiti podložku lepiacou páskou/ Olepit podložku lepící páskou/ Пластика с лентата за запечатване/ Katta plaat isoleerikleerplindiga/ Κολληστε το πλακίδιο με κολλητική ταινία/ Aoperiti placa cu o bandă adezivă/ Prelepiti ploščo/ Peitã mikrotitrauslevy oheisella teipillä
MEASURE		Measure plate within 30 min at 450 nm (Referencefilter ≥590nm)/Ausmessung innerhalb von 30 min bei 450 nm (Referenzfilter ≥ 590 nm)./ Measure lábsorbance en l'éspace de 30 min à 450 nm avec ≥590nm longueur d'onde pour référence/Misurazione entro 30 min. a 450 nm (filtro di riferimento ≥ 590 nm)./ Medición de la placa dentro de los siguientes 30 min a 450 nm (filtro de referencia ≥ 590nm)/ Medir a placa dentro de 30 min a 450 nm (Filtro de referência ≥590nm)/ Binnen 30 minuten bij 450 nm meten (referentiefilter ≥ 590 nm)./ Mål plate i løbet af 30 min ved 450 nm (referencefilter ≥590nm)/ Mät inom 30 min vid 450 nm (referensfilter ≥ 590 nm)./ Pomiar w ciągu 30 min przy 450 nm (filtr odniesienia ≥ 590 nm)./ Ki mérése 30 percen belül 450 nm-nél (referenciaszűrő ≥ 590 nm)./ Merat' 30 minut pri 450 nm/Měřit 30 minut při 450 nm/ Отчитане в рамките на 30 min при 450 nm (референтен филтър ≥ 590 nm)./ Mõõtmise 30 min jooksul 450 nm korral (võrdlusfilter ≥ 590 nm). Μέτρηση εντός 30 min στα 450 nm (φίλτρο αναφοράς ≥ 590 nm)./ Măsurare în decurs de 30 min la 450 nm (filtru de referință ≥ 590 nm)./ Izmerite ploščico v 30 min pri 450 nm (referenčni filter ≥590nm) / Mittaa 30 minuutin aikana 450 nm:ssä (referenssi suodatin ≥ 590 nm)
Literatur		Literature/ Literatur/ Bibliographie/ Letterario/ Bibliografía/ Literatura documentação/ literatuur/ Litteratur/ litteratur/ Literatura/ Irodalom/ Literatúra/ Literatura/ Литература// Kirjandus/ Βιβλιογραφία/ Bibliografie/ literatura/ Lähdeluettelo
International Test description		International test description/ internationale Testanleitung/ description internationale de test/ Istruzioni per il test internazionali/ Descripción de ensayo internacional/ Descrição internacional do teste/ internationale testbeschrijving/ internationell testbeskrivning/ Opis testu międzynarodowego/ nemzetközi teszt-útmutató/ Medzinárodný návod k testu/ Mezinárodní návod k testu/ rahvusvaheline katse kirjeldus/ Διεθνείς οδηγίες για εργαστηριακές δοκιμές/ instructiuni internaționale pentru testare/ mednarodna navodila za preizkus/ Kansainvälinen käyttöohje
End		in all required wells/ in allen benötigten Vertiefungen/ dans tous les godets requis/ in tutti i pozzetti richiesti/ en todos los pozos requeridos/ em todos os tubos necessários/ in alle nodige putjes/ i alle nødvendige brønde/ i alla nödvändiga brunnar/ we wszystkich potrzebnych wgłębieniach/ minden szükséges forrásban/ vo všetkých potrebných miestach/ ve všech potřebných místech/ във всички необходими ямки/ kõigis vajalikes süvendites/ σε όλες τις απαραίτητες κοιλότητες/ în toate cavitățile necesare/ v vseh zahtevanih vdolbinah /kaikkiin tarvittaviin mikrotitrauslevyn syvennyksiin

For Research Use Only.

Not for use in diagnostic procedures.

CAUTION: Not for human or animal therapeutic or diagnostic use.

For in vitro use only.

For professional use only.

Read entire protocol before use!

TABLE OF CONTENTS

INTENDED USE	5
INTRODUCTION	5
REAGENTS PROVIDED	6
MATERIALS REQUIRED BUT NOT PROVIDED.....	6
WARNINGS AND PRECAUTIONS	7
METHOD	8
SPECIMEN	8
STORAGE OF THE SAMPLES.....	8
SAMPLE PREPARATION.....	8
TECHNICAL NOTES.....	8
ASSAY PROCEDURE.....	10
CALCULATION OF RESULTS	11
PERFORMANCE CHARACTERISTICS.....	12
EXEMPLARY VALUES	16
LITERATURE	18
SUMMARY OF THE ASSAY PROCEDURE	19
INTERNATIONALE ASSAY DESCRIPTION.....	20

Mediagnost **Adiponectin ELISA E09**

- **For research and professional use only!**
- is suited for Adiponectin determination in **Serum** and **Plasma** samples,
- is extremely sensitive (less than **0.3 ng/ml**)
- is **fast**: incubation time a total of 1 hour and 45 minutes
- Single Standards with **2, 10, 30, 70, 100 ng/ml** human Adiponectin are provided in the Kit
- 2 Control Sera are provided for quality control purposes according GLP
- is calibrated with native **Adiponectin**
- uses **high affinity monoclonal antibodies** against human Adiponectin
- Microtiter plates are separately breakapart

Intended Use

Measurement of human Adiponectin in human serum and plasma samples for research use!

Introduction

Adiponectin is a 30kDa protein which percentage in serum proteins is 0.01%. It is mainly synthesized by Adipocytes, but also muscle cells and hepatocytes have the ability to synthesize Adiponectin. Until now, IGF-I is the only known natural inductor of the synthesis. It consists of a Collagen-like N-terminal and a globular C-terminal domain [1]. In vivo Adiponectin appears with different oligomers. Beside the trimer and dimer also high molecular multimers exist [1-3]. Up to now two different receptors are known, both receptors are ubiquitarily expressed, though the distribution in the tissues varies. The Adiponectin Receptor 1 (AdipoR1) is especially in muscle- and AdipoR2 in liver tissue synthesized [4].

The significance for the human organism is not clear until now. First studies show, that adiponectin correlates negatively with BMI and thus it could have relevance for the energy metabolism for example through the regulation of fatty acid oxidation. Beside the correlation with BMI, Adiponectin level is associated with the Insulin-Resistance [5-7] and so also linked with Type II Diabetes. Adiponectin is associated also with glucose- und lipometabolism [8, 9].

Furthermore it is involved in inflammatory processes [11-15] and therewith it is of importance for appearance of arteriosclerosis [4, 5, 16] and coronaritis [17, 18], thus the determination of Adiponectin level in plasma could serve to estimate the risk of coronary disease [19, 20]. Beside this Adiponectin influences further physiological processes as for example the angiogenesis [21, 22].

Reagents Provided

1)	MTP	Microtiter plate , ready for use: Microtiter plate with 96 wells, dived up in 12 stripes à 8 wells (separately breakapart), coated with anti-human Adiponectin antibody.
2)	STD	Standards A-E , lyophilised, contain native Adiponectin. Standard values are between 2-100 ng/ml (2, 10, 30, 70 and 100 ng/ml) Adiponectin and have to be reconstituted in 750 µl (each) Dilution Buffer VP . 100 µl per well are used in the assay. If the standards are required for more than one assay process we recommend to store the reconstituted Standards frozen at -20°C. Attention: Standards should be thawed only once – where required please store aliquoted in adequate volumes.
3)	BUF	Dilution buffer VP , 125 ml , ready for use, please use for the reconstitution of Standards A-E, Control Sera KS1 & KS2 and for the serum dilution.
4)	Control	Control Serum KS1 & KS2 , each 500 µl lyophilised: Contains human Serum and has to be reconstituted in 500 µl Dilution Buffer VP . The Adiponectin target value concentration and the respective range are given on the certificate. The dilution of the Control Sera KS 1&2 should be according to the dilution of the respective samples, the target value concentration should be obtained by multiplication with the respective dilution factor .
5)	AbCONJ	Antibody-POD-Conjugate AK , 12 ml , ready for use , contains a mixture of biotinylated anti-Adiponectin antibody and HRP (Horseradish Peroxidase)-labelled Streptavidin. Use 100 µl/well in the assay.
6)	WASHBUF 20x	Washing Buffer (WP) , 50 ml , 20 X concentrated solution. Washing Buffer (WP) has to be diluted 1:20 with distilled or demineralised water before use (e.g. add the complete contents of the flask (50 ml) into a graduated flask and fill up with A.dest. to 1000 ml). Attention: After dilution the Washing Buffer is only 4 weeks stable, dilute only according to requirements.
7)	SUBST	Substrate (S) , 12 ml , ready for use, horseradish-peroxidase-(HRP)-substrate, stabilised H ₂ O ₂ Tetramethylbencidine.
8)	H ₂ SO ₄	Stopping Solution (SL) , 12 ml , ready for use, 0.2 M sulphuric acid, Caution acid!
9)		Sealing tape for covering of the microtiter plate, 2 x, adhesive.

Materials Required but not provided

Precision pipettes Micropipettes and multichannel pipettes with disposable plastic tips

Distilled or deionized water for dilution of the Washing Buffer (WP)

Vortex-mixer

Microtiter plate shaker (350 rpm)

Microtiter plate washer (recommended)

Micro plate reader ("ELISA-Reader") with filter for 450 and ≥590 nm

Polyethylene PE/Polypropylene PP tubes for dilution of samples

WARNINGS AND PRECAUTIONS

For Research Use Only. For Professional use only.

Not for internal use in humans and animals. Follow strictly the test protocol. Use the valid version of the package insert provided with the kit. Be sure that everything has been understood. Mediagnost will not be held responsible for any loss or damage (except as required by statute) howsoever caused, arising out of noncompliance with the instructions provided. A Material Safety Data sheet is available on request.

Do not use obviously damaged or microbial contaminated or spilled material.

Use separate pipette tips for each sample, control and reagent to avoid cross contamination. Use reservoirs only for single reagents. This especially applies to the substrate reservoirs. Using a reservoir for dispensing a substrate solution that had previously been used for the conjugate solution may turn solution colored. Do not pour reagents back into vials as reagent contamination may occur. Mix the contents of the microplate wells thoroughly to ensure good test results. Do not reuse microwells. Do not let wells dry during assay; add reagents immediately after completing the rinsing steps.

Caution: This kit contains material of human and/or animal origin. Therefore all components and specimens should be treated as potentially infectious.

Appropriate precautions and good laboratory practices must be used in the storage, handling and disposal of the kit reagents. The disposal of the kit components must be made according to the local regulations.

Human Serum

Following components contain human serum: **Control Sera KS1 and KS2, and Standards A.E**

Source human serum for the control sera provided in this kit was tested and found non-reactive for Hepatitis-B surface antigen (HBsAg), Hepatitis C virus (HCV), and Human Immunodeficiency Virus 1 and 2 (HIV). No known methods can offer total security of the absence of infectious agents; therefore all components and specimens should be treated as potentially infectious.

Reagents A-E, AK, VP, WP

Contain as preservative a mixture of **5-chloro-2-methyl-4-isothiazolin-3-one and 2-methyl-4-isothiazolin-3-one** (<0.015%)

H317	May cause an allergic skin reaction.
P280	Wear protective gloves/ protective clothing/ eye protection/ face protection.
P272	Contaminated work clothing should not be allowed out of the workplace.
P261	Avoid breathing dust/ fume/ gas/ mist/ vapours/ spray.
P333+P313	If skin irritation or rash occurs: Get medical advice/ attention.
P302+P352	IF ON SKIN: Wash with plenty of soap and water.
P501	Dispose of contents/ container in accordance with local/ regional/ national/ international regulations.

Substrate Solution (S)

The TMB-Substrate (S) contains 3,3',5,5' Tetramethylbencidine (<0.05%)

H315	Causes skin irritation.
H319	Causes serious eye irritation.
H335	May cause respiratory irritation.
P261	Avoid breathing dust/ fume/ gas/ mist/ vapours/ spray.
P305+P351+	IF IN EYES: Rinse cautiously with water for several minutes.
P338	Remove contact lenses, if present and easy to do. Continue rinsing.

Stopping Solution (SL)

The Stopping solution contains 0.2 M acid sulphur acid (H₂SO₄)

H290	May be corrosive to metals.
H314	Causes severe skin burns and eye damage.
P280	Wear protective gloves/ protective clothing/ eye protection/ face protection.
P301+P330+	IF SWALLOWED: rinse mouth.
P331	Do NOT induce vomiting.
P305+P351+	IF IN EYES: Rinse cautiously with water for several minutes.
P338	Remove contact lenses, if present and easy to do. Continue rinsing.
P309+P310	IF exposed or if you feel unwell: Immediately call a POISON CENTER or doctor/physician.

General first aid procedures:

Skin contact: Wash affected area rinse immediately with plenty of water at least 15 minutes. Remove contaminated cloths and shoes.

Eye contact: In case of contact with eyes, rinse immediately with plenty of water at least 15 minutes. In order to assure an effectual rinsing spread the eyelids.

Ingestion: After swallowing the product, if the affected person is conscious, rinse out the mouth with plenty of water: seek medical advice immediately.

Method

The Mediagnost ELISA for Adiponectin E09 is a so-called Sandwich-Assay using two specific and high affine antibodies. The Adiponectin in the samples binds to the first antibody coated on the microtiter plate. In the following step the second specific anti-Adiponectin-Antibody binds in turn to the immobilised Adiponectin. The second antibody is biotinylated and will be applied in a mixture with a Streptavidin-Peroxidase-Enzyme Conjugate. In the closing substrate reaction the turn of the colour will be catalysed depending on the Adiponectin-level of the samples.

Specimen

Serum and heparin plasma levels are comparable. In EDTA- and Citrate Plasma-samples levels were found approx. 18% lower, because of the relatively high amount of anticoagulant. Adiponectin can also be measured in urine, breast milk and cellculture media by this testsystem.

The blood sample for serum preparation should be gained according to standardized venipuncture procedure. The samples should be stored without anticoagulation reagents. Hemolytic reactions have to be avoided. The blood has to be allowed to clot and after complete clotting, serum is separated by centrifugation.

Storage of the samples

Storage at RT max. 2 days

Storage at -20°C max. 2 years

More than three freeze/thaw cycles are not possible.

Sample Preparation

Samples have to be diluted in Dilution Buffer (VP). The excellent linearity of this testsystem allows sample dilution of 1:200 to 1:1600.

We recommend a standard dilution of 1:310.

Suggestion for dilution protocol:

Dilute for example 300 μl Dilution Buffer VP in PE-/PP-Tubes (application of a multi-stepper is recommended in larger series), add 10 μl Serum- or Plasma (dilution: 1:31). Add **900 μl** Dilution Buffer **VP** in an other PE-/PP-tube and **100 μl** of the thoroughly mixed first dilution. After mixing, use **2x100 μl** from this **1:310** diluted sample in the assay.

If you have the necessary technical equipment a one-step dilution of 1:301 is possible: Add 5 μL to 1.5 mL dilution buffer VP.

Technical Notes

The assay has to be conducted strictly according the test protocol herein.

Reagents with different lot numbers cannot be mixed. The microtiterplate and reagents are stable until the indicated expiry if stored unopened and protected from sunlight at $2 - 8^{\circ}\text{C}$.

The shelf life of the components after initial opening is warranted for 4 weeks, store the unused strips and microtiter wells airtight together with the desiccant at $2-8^{\circ}\text{C}$ in the clip-lock bag, use in the frame provided.

Bring all reagents to room temperature ($20 - 25^{\circ}\text{C}$) before use. Possible precipitations in the buffers have to be resolved before usage by mixing and / or warming.

Use separate pipette tips for each sample, control and reagent to avoid cross contamination. Use reservoirs only for single reagents. This especially applies to the substrate reservoirs. Using a reservoir for dispensing a substrate solution that had previously been used for the conjugate solution may turn solution colored. Do not pour reagents back into vials as reagent contamination may occur. Mix the contents of the microplate wells thoroughly to ensure good test results. Do not reuse microwells. Do not let wells dry during assay; add reagents immediately after completing the rinsing steps.

Incubation at room temperature means: 20-25°C

The incubation steps should be performed at mean rotation frequency of a particularly suitable microtiter plate shaker. We are recommending 350 rpm. Due to certain technical differences deviations may occur, in case the rotation frequency must become adjusted. Insufficient shaking may lead to inadequate mixing of the solutions and thereby to low optical densities, high variations and/or false values, excessive shaking may result in high optical densities and/or false values.

Proper washing is of basic importance for a secure, reliable and precise performance of the test. Incomplete washing is common and will adversely affect the test outcome. Possible consequences may be uncontrolled unspecific variations of measured optical densities, potentially leading to false results calculations of the examined samples. Effects like high background values or high variations may indicate washing problems.

All washing must be performed with the provided washing buffer diluted to usage concentration. Washing volume per washing cycle and well must be 300 µl at least.

The danger of handling with potentially infectious material must be taken into account.

When using an automatic microtiter plate washer, the respective instructions for use must be carefully followed. Device adjustments, e.g. for plate geometry and the provided washing parameters, must be performed. Dispensing and aspirating manifold must not scratch the inside well surface. Provisions must be made that the remaining fluid volume of every aspiration step is minimized. Following the last aspiration step of each washing cycle, this could be controlled, and possible remaining fluid could then be removed, by inverting the plate and repeatedly tapping it dry on non fuzzy absorbent tissue.

Manual washing is an adequate alternative option. Washing Buffer may be dispensed via a multistepper device, a multichannel pipette, or a squirt bottle. The fluid may be removed by dynamically swinging out the microtiter plate over a basin. If aspirating devices are used, care has to be taken that the inside well surface is not scratched. Subsequent to every single washing step, the remaining fluid should be removed by inverting the plate and repeatedly tapping it dry on non fuzzy absorbent tissue.

Standards and Controls

For the reconstitution of the lyophilised components (Standards A - E and Control Sera KS1 &KS2) the kit Dilution Buffer VP has to be used. It is recommended to keep reconstituted reagents at room temperature for 15 minutes and then to mix them thoroughly but gently (no foam!) with a Vortex mixer.

The reconstituted standard and controls can be stored for 1 month at -20°C. Repeated freeze/thaw cycles have to be avoided.

Washing Buffer

The required volume of washing buffer is prepared by 1:20 dilution of the provided 20fold concentrate with deionised water. The diluted Washing Buffer is stable for 4 weeks at 2-8°C. It has to be at room temperature for usage!

Assay Procedure

All determinations (Standards, Control Sera KS1 & KS2 and samples) should be assayed in duplicate. For optimal results, accurate pipetting and adherence to the protocol are recommended.

When performing the assay, the Standards, Control Sera and the samples should be pipetted as fast as possible (e.g., <15 minutes).

All incubations have to be conducted at room temperature (20-25°C)

To avoid distortions due to differences in incubation times, Antibody-POD-Conjugate **AK** as well as the following Substrate Solution **S** should be added to the plate in the same order and in the same time interval as the samples. Stop Solution **SL** should be added to the plate in the same order as the Substrate Solution.

- 1) Add **100 µl Dilution Buffer VP** in the first wells. Subsequently add **100 µl Standard** or **100 µl** of diluted **Control Sera** or **diluted samples**.
- 2) Cover the wells with sealing tape and incubate the plate for **1 hour at room temperature** (shake at 350 rpm)
- 3) After incubation aspirate the contents of the wells and wash the wells **3 times 300 µl Washing Buffer WP** / well. The washing buffer should incubate for at least for 15 seconds/cycle.
- 4) Following the last washing step pipette **100 µl** of the **Antibody-POD-Conjugate AK** in each well.
- 5) Cover the wells with sealing tape and incubate the plate for **30 Minutes at room temperature** (shake at 350 rpm).
- 6) After incubation wash the wells **3 times** with **Washing Buffer WP** as described in step 3
- 7) Pipette **100 µl of the TMB Substrate Solution S** in each well.
- 8) Incubate the plate for **15 minutes in the dark at room temperature (20 - 25°C)**.
- 9) Stop the reaction by adding **100 µl of Stopping Solution SL**.
- 10) Measure the colour reaction within 30 minutes at 450nm (reference filter ≥590 nm).

Calculation of Results

Establishing the Standard Curve

For the evaluation of the assay it is preconditioned that the absorbance values of the blank should be below 0.25, these of standard E should be above 1.0.

Samples, which yield higher absorbance values than Standard E are beyond the standard curve, for reliable determinations these samples should be tested anew with a higher dilution.

Standards are provided in the following concentrations (use the concentration unit as preferred):

Standard	A	B	C	D	E
ng/ml	2	10	30	70	100
µg/ml	0.002	0.01	0.03	0.07	0.1

- 1) Calculate the **mean absorbance** value for the blank from the duplicated determination (well A1/A2).
- 2) Subtract the mean absorbance of the blank from the mean absorbances of all other values.
- 3) Plot the standard concentrations on the x-axis versus the mean value of the absorbance of the standards on the y-axis.
- 4) Recommendation: Calculation of the standard curve should be done by using a computer program. **A higher-grade polynomial, or four parametric logistic (4-PL) curve fit or non-linear regression** are usually suitable for the evaluation (as might be spline or point-to-point alignment in individual cases).
- 5) The Adiponectin concentration of the diluted sample or the diluted control sera KS1&KS2 in ng/ml (or µg/ml according the chosen unit for the standards) is calculated in this way, the Adiponectin concentration of the **undiluted sample** and of KS1 & KS2 is calculated **by multiplication** with the respective dilution factor.

The exemplary shown standard curve in Fig.1 **cannot** be used for calculation of your test results. You have to establish a standard curve for each test you conduct!

Exemplary calculation of the adiponectin concentration of a 1:310 diluted sample:

Measured extinction of your sample 0.39
Measured extinction of the blank 0.04

Your measurement programm will calculate the adiponectin concentration of the diluted sample automatically by using the difference of sample and blank for the calculation. You only have to determine the most suitable curve fit (here: polynomial 3rd degree).

In this exemplary case the following equation is solved by the programm to calculate the adiponectin concentration in the sample:

$$0.35 = 5 \times 10^{-7} x^3 - 0.0002x^2 + 0.0346x - 0.0166$$
$$11.13 = x$$

If the dilution factor (**1:310**) is taken into account the adiponectin concentration of the undiluted sample is

$$11.13 \times 310 = 3450.3 \text{ ng/mL} = 3.45 \text{ µg/mL}$$

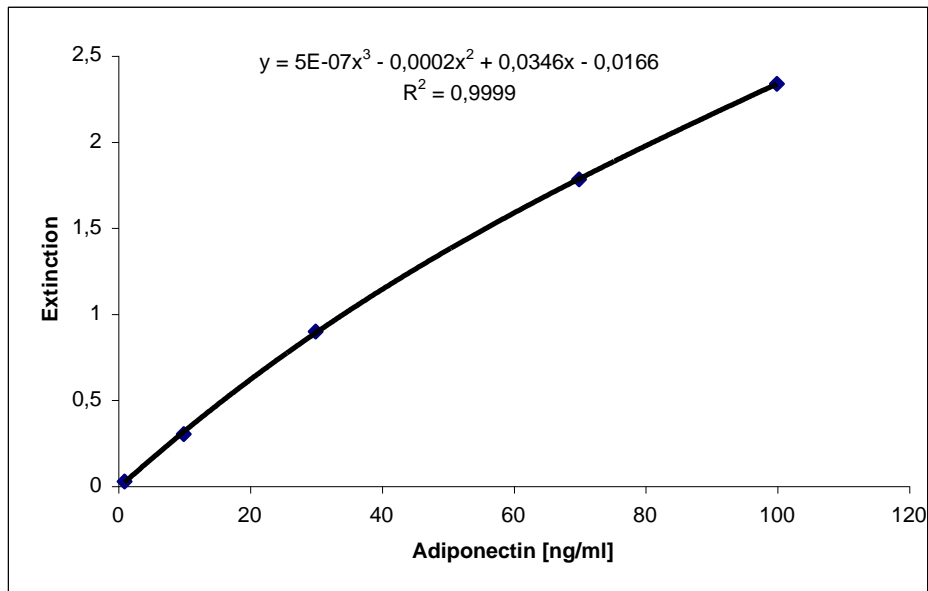


Figure 1 Exemplary Standard Curve with a polynomial 3rd degree as curve fit.

Performance Characteristics

Standards

The Standards of the ELISA E09 are prepared from **native Adiponectin** (Human Serum) in concentrations of 2, 10, 30, 70 and 100 ng/ml. The native Adiponectin was quantified with a recombinant protein and with a commercialised radio immunoassay (Linco Corp.) for Adiponectin.

Sensitivity

The analytical sensitivity of the ELISA E09 was measured by the variability of the signal of the null standard. Based on the twofold standard deviation of the blank the mean analytical sensitivity is < 0.27 ng/mL (Range 0.094 to 0.59 ng/mL).

Specificity

Adiponectin exists in different oligomeric forms: the high, medium and low molecular weight form. Different numbers of the adiponectin monomer aggregate specifically to form a complex. In Figure 2a the five different forms of human adiponectin are shown schematically. In parallel the results of a size-exclusion chromatography of human serum measured with the Mediagnost E09 Adiponectin ELISA are shown.

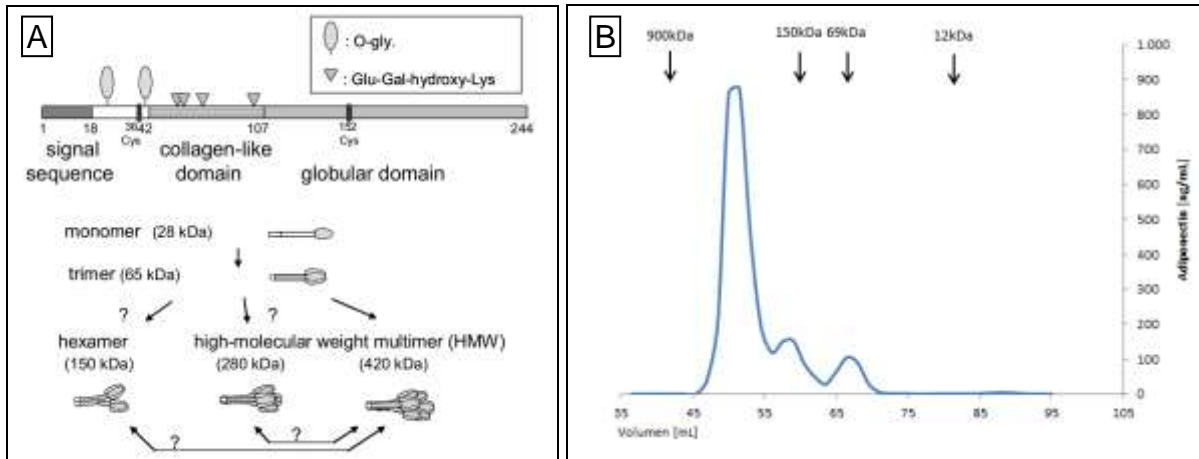


Figure 2 Adiponectin Structure (A) proteinaceous structure of human adiponectin including posttranslational modifications and multimeric forms (taken from Nakano et al¹). (B) Results of a serum separation by size-exclusion chromatography. The sample was fractionated and the adiponectin content of each fraction was measured by Mediagnost E09 Adiponectin, for comparison the corresponding size in kDa is indicated by black arrows.

The results shown in Figure 2b clearly demonstrate that the Mediagnost E09 Adiponectin ELISA detects all forms of Adiponectin present in human serum: the trimer at 65 kDa, the hexamer at 150 kDa and the high molecular weight forms of >280 kDa. The Mediagnost E09 Adiponectin ELISA therefore measures total Adiponectin.

Serum of the cited species was diluted and used as sample in this assay system. No cross reactivity was detected:

Horse, Cow, Chicken, Rabbit, Dog, Guinea pig, Sheep, Mouse, Goat, Donkey, Rat, Cat
Whether this obvious non-reactivity is species specific should be assessed individually by each customer.

Interference

Interference of bilirubin and triglycerides was tested by adding different amounts of these substances to human serum containing adiponectin. For comparison the same amount of buffer without any substance was also added to the serum. Table 1 demonstrates that neither bilirubin nor triglycerides or haemoglobin exert any influence on the measurement of adiponectin in human serum.

Table 1 Interference. Serum samples were enriched with different amount of triglycerides, bilirubin or haemoglobin. The relative amount of adiponectin measured in comparison with native serum is shown here [%].

Triglyceride 100 mg/mL	Bilirubin 100 µg/mL	Hemoglobin 1 mg/mL	Hemoglobin 5 mg/mL
94	96	90	109
90	93	97	--
95	94	93	--

¹ Nakano Y, Tahima S, Yoshimi A, Akiyama H, Tsushima M, Tanioka T, Negoro T, Tomita M, Tobe T: A novel enzyme-linked immunosorbent assay specific for high-molecular-weight adiponectin. J Lipid Res 2006, 47(7):1572-1582

Recovery

Trueness and traceability of the Mediagnosot Adiponectin ELISA E09 was evaluated by recovery of recombinant adiponectin in human serum. The recovery of recombinant Adiponectin yielded in a serum matrix on average 110%.

Table 2 Recovery of recombinant human Adiponectin in Serum. Recombinant Adiponectin was added in different amounts to human serum. The Adiponectin content of the so enriched samples was measured and recovery in comparison to enriched buffer calculated.

R&D recombinant Adiponectin lot 1022911	VP	Serum	Recovery
ng/mL	ng/mL	ng/mL	%
0	0	0.00778	---
75	87.92	95.71	109
37.5	51.84	59.98	116
18.75	26.55	26.7	101
9.375	13.35	15.15	113

Precision

Intra-Assay Variance & Accuracy

Intra assay variance and accuracy is exemplarily shown with two samples. The adiponectin concentration of these samples was repeatedly measured in one assay.

Table 3 Intra-Assay Variation. Recombinant adiponectin was diluted in dilution buffer and the adiponectin concentration of the dilution was measured repeatedly within one assay.

	Determinations [n]	Mean value [µg/L]	Standard deviation [µg/L]	VC [%]	Target Value [µg/L]
Sample 1	8	7.108	0.22	3.14	6
Sample 2	8	107.96	3.97	3.67	100

In both samples the variability is less than 5% and the deviation from the target value is <20%.

Inter-Assay Variance

Serum samples were repeatedly measured in independent assays of different lots. On average the coefficient of variation was 7.5% (SD 1.6). The results of 5 samples are shown in table 4.

Table 4 Inter-Assay Variation

	Sample 1	Sample 2	Sample 3	Sample 4	Sample 5
Mean [µg/mL]	4.72	5.25	8.36	5.45	22.29
CV [%]	8.16	8.14	6.93	8.05	7.30
n	99	68	62	174	62

Linearity

Linearity of sample dilution was tested by serial dilution (1:100 – 1:4000) of human serum samples and recalculation of the adiponectin content in comparison to the mean adiponectin concentration of all dilutions (Table 5), no diluted sample showed a deviation of >30%.

Table 5 Linearity. Human serum sample were diluted in VP and adiponectin content was recalculated. Measurements results are shown in [mg/L]. No deviation of the mean >30% was detected.

µg/mL	Mean	1:100	1:200	1:400	1:800	1:1000	1:2000	1:4000
Sample 1	5.76	6.53	6.331	5.764	5.49	6.067	6.114	4.056
Sample 2	11.53	10.93	12.107	11.395	11.454	11.567	12.884	10.362
Sample 3	12.07	13.57	12.853	12.03	11.974	11.338	11.548	11.169
Sample 4	4.89	4.659	4.886	4.384	4.425	5.851	5.13	n/a

Additionally dilutions of 1:1500 to 1:96000 were evaluated with two samples. In Figure 3 the results are shown and demonstrate that in the tested samples no effect of dilution could be detected on measured adiponectin concentrations. The deviation of the target concentration of each dilution was less than 30% in all dilutions.

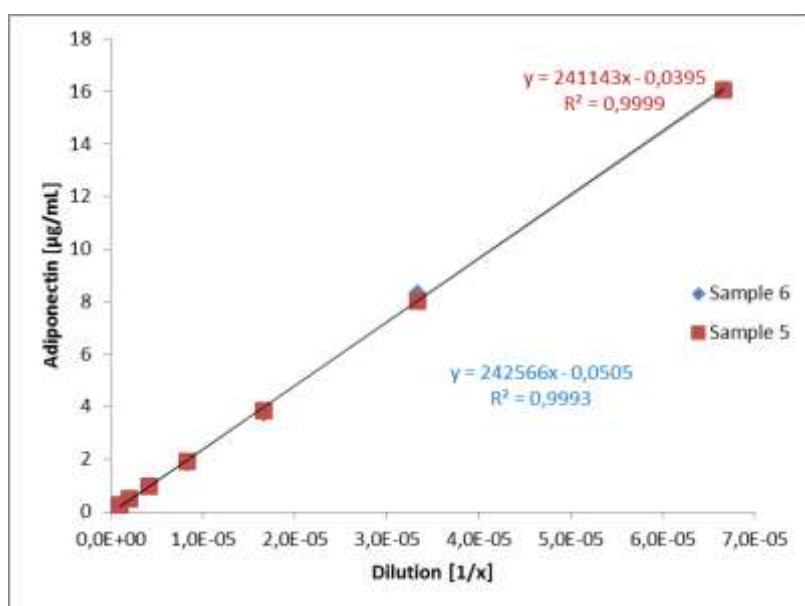


Figure 3 Linearity. Adiponectin concentration was measured in two human samples diluted 1:1500 to 1:96000.

Exemplary Values

The exemplary values for serum adiponectin, were measured with the Mediagnost ELISA E09 in healthy donors and analysed by Prof. Dr. J. Kratzsch, Department of Laboratory Medicine, University Hospital Leipzig, are given below (Tab. 6).

These data show significant correlation between Adiponectin-Serum values and age as well as gender of the probands, in turn the correlation between the respective BMI seems to be less significant. In the samples of neonatal cord blood very high values were found.

Table 6a Exemplary values for **adults, gender specific** mean as well as median, 5. and 95. percentiles are given.

Sex	number	Mean [µg/ml]	Median [µg/ml]	Standard Deviation	5 th Percentile [µg/ml]	95 th Percentile [µg/ml]
Female	101	10.2	9.1	4.6	4.0	19.4
Male	125	6.8	6.1	4.1	2.0	13.9
total	226	8.3	7.5	4.6	2.4	19.3

Table 6b Exemplary values for **children, gender specific** mean as well as median, 5. and 95. percentiles are given.

Sex	number	Mean [µg/ml]	Median [µg/ml]	Standard Deviation	5 th Percentile [µg/ml]	95 th Percentile [µg/ml]
Female	131	8.71	8.18	4.32	3.05	15.6
Male	134	8.97	8.12	5.13	3.36	18.6
total	265	8.84	8.18	4.74	3.33	16.5

Table 6c Exemplary values for Adiponectin, **age specific** mean as well as median, 5. and 95. percentiles are given.

Age (in years)	number	Mean [µg/ml]	Median [µg/ml]	5 th Percentile [µg/ml]	95 th Percentile [µg/ml]
< 7.99	46	12.82	11.7	2.33	26.5
8 – 9.99	40	8	8.09	3.96	14.9
10-11.99	55	8.02	7.14	3.36	13.8
12 – 13.99	26	8.21	7.54	4.5	13.2
14 – 15.99	59	8.12	8.09	3.67	13.7
16 – 19.99	41	7.97	7.79	2.74	13.3
all	267	8.88	8.18	3.33	16.7
20 – 29.99	47	6.72	6.38	2.5	12.25
30 – 39.99	38	7.38	6.69	1.98	19.29
40 – 49.99	55	8.42	8.20	2.41	17.85
50 – 59.99	47	9.61	8.55	2.15	19.85
> 60	32	9.52	8.57	3.00	21.10
all	226	8.33	7.5	2.41	19.29

Table 6d Exemplary values for Adiponectin, **age** as well as **gender specific** mean and median, BMI and 25. and 75. percentiles are given.

Female			Adiponectin (µg/ml):			
Age (Years):	n:	BMI: AV ± SD	AV ± SD::	Median :	Percentile: 25.- 75.	Min. – Max.:
Newborn Cord blood	19		29.80 ± 12.49	26.1	19.5-35.2	16.9-61.4
< 3.99	9	15.73 ± 0.79	14.43 ± 7.76	11.2	8.2-21.8	2.3-26.7
4.0-7.99	11	16.01 ± 1.94	8.46 ± 4.73	9.3	2.9-12.1	1.4-15.6
8.0-9.99	22	17.58 ± 3.84	7.92 ± 3.00	8.2	5.2-10.0	3.6-15.1
10.0-11.99	33	17.83 ± 1.86	7.66 ± 4.59	6.6	5.0-8.8	3.1-20.9
12.0-13.99	11	19.85 ± 2.31	8.22 ± 5.64	7.5	6.5-9.2	4.9-13.2
14.0-15.99	27	19.91 ± 1.72	8.83 ± 9.25	8.9	5.2-11.8	2.6-17.7
16.0-19.99	18	21.64 ± 2.64	9.00 ± 3.22	8.7	6.9-11.2	2.7-14.0
20.0-29.99	24	23.12 ± 5.01	7.39 ± 3.35	7.3	5.7-9.0	3.4-17.8
30.0-39.99	17	23.20 ± 2.86	9.19 ± 3.89	8.6	7.2-10.4	3.6-19.3
40.0-49.99	26	24.50 ± 4.11	9.93 ± 3.59	9.5	7.5-11.6	4.4-19.6
50.0-59.99	21	24.61 ± 3.31	11.5 ± 5.49	10.0	8.0-15.9	2.0-23.1
>60.0	8	24.63 ± 1.89	15.6 ± 4.64	15.3	11.4-18.2	11.2-24.1

n= Number of Proband **AV**=Average Value, **BMI**=Body Mass Index (kg/m²) **SD**=Standard Deviation

Male			Adiponectin (µg/ml):			
Age (Years):	n:	BMI: AV ± SD	AV ± SD:	Median :	Percentile: 25.- 75.	Min. – Max.:
Newborn Cord blood	10		27.80 ± 7.68	26.7	22.2-31.0	15.6-40.6
< 3.99	14	16.17 ± 1.81	16.57 ± 6.55	14.3	11.6-21.2	5.8-40.3
4.0-7.99	12	15.69 ± 1.05	11.24 ± 5.43	9.7	8.9-15.9	3.5-18.6
8.0-9.99	18	16.45 ± 1.76	8.11 ± 2.93	7.6	6.2-9.1	5.00-15.4
10.0-11.99	21	18.34 ± 2.18	8.43 ± 3.91	7.8	5.2-10.9	3.4-20.2
12.0-13.99	14	18.61 ± 2.11	7.59 ± 2.86	7.1	6.0-10.3	2.4-12.2
14.0-15.99	32	19.86 ± 2.00	7.53 ± 2.52	7.4	5.1-9.3	3.8-15.4
16.0-19.99	23	22.03 ± 2.42	7.16 ± 3.53	6.9	4.2-9.6	2.0-13.9
20.0-29.99	23	23.43 ± 2.48	5.44 ± 2.29	5.8	4.0-6.9	1.3-10.3
30.0-39.99	21	23.33 ± 2.72	5.92 ± 4.60	4.4	2.7-6.7	1.9-20.6
40.0-49.99	22	23.79 ± 2.41	6.13 ± 2.92	5.5	3.8-8.3	2.1-11.6
50.0-59.99	23	26.68 ± 2.77	7.45 ± 4.50	6.7	5.0-8.8	1.4-19.6
>60.0	24	25.72 ± 2.12	7.48 ± 3.92	7.6	4.6-9.2	3.0-21.1

n= Number of Proband **AV**=Average Value, **BMI**=Body Mass Index (kg/m²) **SD**=Standard Deviation


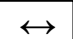

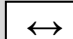


Literature

1. Nakano, Y., et al., *Isolation and characterization of GBP28, a novel gelatin-binding protein purified from human plasma*. J Biochem (Tokyo), 1996. **120**(4): p. 803-12.
2. Pajvani, U.B., et al., *Structure-function studies of the adipocyte-secreted hormone Acrp30/adiponectin. Implications for metabolic regulation and bioactivity*. J Biol Chem, 2003. **278**(11): p. 9073-85.
3. Tsao, T.S., et al., *Role of disulfide bonds in Acrp30/adiponectin structure and signaling specificity. Different oligomers activate different signal transduction pathways*. J Biol Chem, 2003. **278**(50): p. 50810-7.
4. Shimada, K., T. Miyazaki, and H. Daida, *Adiponectin and atherosclerotic disease*. Clin Chim Acta, 2004. **344**(1-2): p. 1-12.
5. Higashiura, K., et al., *Correlations of adiponectin level with insulin resistance and atherosclerosis in Japanese male populations*. Clin Endocrinol (Oxf), 2004. **61**(6): p. 753-9.
6. Spranger, J., et al., *Adiponectin is independently associated with insulin sensitivity in women with polycystic ovary syndrome*. Clin Endocrinol (Oxf), 2004. **61**(6): p. 738-46.
7. Zoico, E., et al., *Adipocytokines, fat distribution, and insulin resistance in elderly men and women*. J Gerontol A Biol Sci Med Sci, 2004. **59**(9): p. M935-9.
8. Ye, J.M., et al., *PPARalpha /gamma ragaglitazar eliminates fatty liver and enhances insulin action in fat-fed rats in the absence of hepatomegaly*. Am J Physiol Endocrinol Metab, 2003. **284**(3): p. E531-40.
9. Yamauchi, T., et al., *Adiponectin stimulates glucose utilization and fatty-acid oxidation by activating AMP-activated protein kinase*. Nat Med, 2002. **8**(11): p. 1288-95.
10. Blüher (Blueher), M., et al., *Total and high-molecular weight adiponectin in relation to metabolic variables at baseline and in response to an exercise treatment program: comparative evaluation of three assays*. Diabetes Care, 2007. **30**(2): p. 280-5.
11. Delaigle, A.M., et al., *Induction of adiponectin in skeletal muscle by inflammatory cytokines: in vivo and in vitro studies*. Endocrinology, 2004. **145**(12): p. 5589-97.
12. Winzer, C., et al., *Plasma adiponectin, insulin sensitivity, and subclinical inflammation in women with prior gestational diabetes mellitus*. Diabetes Care, 2004. **27**(7): p. 1721-7.
13. Xydakis, A.M., et al., *Adiponectin, inflammation, and the expression of the metabolic syndrome in obese individuals: the impact of rapid weight loss through caloric restriction*. J Clin Endocrinol Metab, 2004. **89**(6): p. 2697-703.
14. Motoshima, H., et al., *Adiponectin suppresses proliferation and superoxide generation and enhances eNOS activity in endothelial cells treated with oxidized LDL*. Biochem Biophys Res Commun, 2004. **315**(2): p. 264-71.
15. Wolf, A.M., et al., *Adiponectin induces the anti-inflammatory cytokines IL-10 and IL-1RA in human leukocytes*. Biochem Biophys Res Commun, 2004. **323**(2): p. 630-5.
16. Okamoto, Y., et al., *Adiponectin reduces atherosclerosis in apolipoprotein E-deficient mice*. Circulation, 2002. **106**(22): p. 2767-70.
17. Schlegel, A., *Adiponectin and risk of coronary heart disease*. Jama, 2004. **292**(1): p. 40; author reply 40.
18. Choi, K.M., et al., *Inflammation, Insulin Resistance, and Glucose Intolerance in Acute Myocardial Infarction Patients without a Previous Diagnosis of Diabetes Mellitus*. J Clin Endocrinol Metab, 2004.
19. Nakamura, Y., et al., *Implications of plasma concentrations of adiponectin in patients with coronary artery disease*. Heart, 2004. **90**(5): p. 528-33.
20. Pischon, T., et al., *Plasma adiponectin levels and risk of myocardial infarction in men*. Jama, 2004. **291**(14): p. 1730-7.
21. Shibata, R., et al., *Adiponectin stimulates angiogenesis in response to tissue ischemia through stimulation of amp-activated protein kinase signaling*. J Biol Chem, 2004. **279**(27): p. 28670-4.
22. Fernandez-Real, J.M., et al., *Adiponectin is associated with vascular function independent of insulin sensitivity*. Diabetes Care, 2004. **27**(3): p. 739-45.

SUMMARY OF THE ASSAY PROCEDURE

Preparation of reagents		Reconstitution:	Dilution
A-E	Standards	in 750 µL Dilution Buffer VP	-
KS1 and KS2	Control Sera	in 500 µL Dilution Buffer VP	1:310 with VP
WP	Washing Buffer	-	1:20 with Aqua dest.
Sample dilution: with Dilution Buffer VP 1:310			
Before assay procedure bring all reagents to room temperature 20-25°C .			
Assay Procedure in Double Determination:			
Pipette	Reagents	Position	
100 µL	Dilution Buffer VP (Blank)	A1/A2	
100 µL	Standard A (2 ng/mL)	B1/B2	
100 µL	Standard B (10 ng/mL)	C1/C2	
100 µL	Standard C (30 ng/mL)	D1/D2	
100 µL	Standard D (70 ng/mL)	E1/E2	
100 µL	Standard E (100 ng/mL)	F1/F2	
100 µL	Control Serum KS1 (1:310 diluted)	G1/G2	
100 µL	Control Serum KS2 (1:310 diluted)	H1/H2	
100 µL	Sample (1:310 diluted)	in the rest of the wells according the requirements	
Cover the wells with the sealing tape.			
Sample Incubation: 1 h at 20-25°C, 350 rpm			
3 x 300 µL	Aspirate the contents of the wells and wash 3 x with 300 µL each Washing Buffer WP/ well	In each well	
100 µL	Antibody-POD-Conjugate AK	In each well	
Cover the wells with the sealing tape.			
Incubation: 30 Minutes at 20-25°C, 350 rpm			
3 x 300 µL	Aspirate the contents of the wells and wash 3 x with 300 µL each Washing Buffer WP/ well	In each well	
100 µL	Substrate Solution S	In each well	
Incubation: 15 Minutes in the Dark at 20-25°C			
100 µL	Stopping Solution SL	In each well	
	Measure the absorbance within 30 min at 450 nm with ≥ 590 nm as reference wavelength.		

Internationale Assay Description

A-E	STD	Rec in 750 µL	BUF VP	-
KS1	Control	Rec in 500 µL	BUF VP	1:310 DILU BUF VP
KS2	Control	Rec in 500 µL	BUF VP	1:310 DILU BUF VP
WP	WASHBUF 20x	-	-	1:20 DILU A. dest.
-	SPE	-	-	1:310 DILU BUF VP
-	°C	20-25 °C		
100 µL	BUF VP			A1/A2
100 µL	STD A	(2 ng/mL)		B1/B2
100 µL	STD B	(10 ng/mL)		C1/C2
100 µL	STD C	(30 ng/mL)		D1/D2
100 µL	STD D	(70 ng/mL)		E1/E2
100 µL	STD E	(100 ng/mL)		F1/F2
100 µL	CONTROL	KS1 1:310	DILU BUF VP	G1/G2
100 µL	CONTROL	KS2 1:310	DILU BUF VP	H1/H2
100 µL	SPE	1:310	DILU BUF VP	
TAPE				
 1 h °C 20-25  350 rpm				
3x 300 µL	3x WASHBUF WP			
100 µL	AbConj AK			
TAPE				
 0.5 h °C 20-25  350 rpm				
3x 300 µL	3x WASHBUF WP			
100 µL	SUBST TMB S			
 0.25 h °C 20-25 				
H ₂ SO ₄ SL				
MEASURE				