

# ORGENTEC Diagnostika GmbH

Carl-Zeiss-Straße 49-51

55129 Mainz - Germany

Phone: +49 (0) 61 31 / 92 58-0

Fax: +49 (0) 61 31 / 92 58-58

Internet: [www.orgentec.com](http://www.orgentec.com)



**For research use only.  
Not for diagnostic procedures.**



## RUO 523 Anti-BPI

### INTENDED USE

Anti-BPI is an ELISA test system for the quantitative measurement of IgG class autoantibodies against bactericidal permeability increasing protein, BPI, in human serum or plasma. This product is intended for research use only. Not for diagnostic procedures.



Distributed By:

**IBL-America, Inc.**

8201 Central Ave NE, Suite P

Minneapolis, MN 55432, USA

[info@ibl-america.com](mailto:info@ibl-america.com)

(888) 523 1246

### SYMBOLS USED ON LABELS



Manufacturer



Catalogue number



Sufficient for 96 determinations



Batch code



Use by



Temperature limitation



Keep away from sunlight



Do not reuse



Date of manufacture



Consult instructions for use

MICROPLATE

Microplate

CALIBRATOR A

Calibrator

CALIBRATOR B

Calibrator

CALIBRATOR C

Calibrator

CALIBRATOR D

Calibrator

CALIBRATOR E

Calibrator

CALIBRATOR F

Calibrator

CONTROL +

Control positive

CONTROL -

Control negative

DILUENT

Sample Buffer P

CONJUGATE

Enzyme Conjugate

TMB

TMB Substrate

STOP

Stop solution

WASH

Wash Buffer

RTU

Ready to use

## PRINCIPLE OF THE TEST

Highly purified BPI is bound to microwells.

The determination is based on an indirect enzyme linked immune reaction with the following steps:

Specific antibodies in the sample bind to the antigen coated on the surface of the reaction wells. After incubation, a washing step removes unbound and unspecifically bound serum or plasma components. Subsequently added enzyme conjugate binds to the immobilized antibody-antigen-complexes. After incubation, a second washing step removes unbound enzyme conjugate. After addition of substrate solution the bound enzyme conjugate hydrolyses the substrate forming a blue coloured product. Addition of an acid stops the reaction generating a yellow end-product. The intensity of the yellow color

correlates with the concentration of the antibody-antigen-complex and can be measured photometrically at 450 nm.

## WARNINGS AND PRECAUTIONS

- All reagents of this kit are intended for professional use only.
- Components containing human serum were tested and found negative for HBsAg, HCV, HIV1 and HIV2 by FDA approved methods. No test can guarantee the absence of HBsAg, HCV, HIV1 or HIV2, and so all human serum based reagents in this kit must be handled as though capable of transmitting infection.
- Bovine serum albumin (BSA) used in components has been tested for BSE and found negative.
- Avoid contact with the substrate TMB (3,3',5,5'-Tetramethyl-benzidine).
- Stop solution contains acid, classification is non-hazardous. Avoid contact with skin.
- Control, sample buffer and wash buffer contain sodium azide 0.09% as preservative. This concentration is classified as non-hazardous.
- Enzyme conjugate contains ProClin 300 0.05% as preservative. This concentration is classified as non-hazardous.

During handling of all reagents, controls and serum samples observe the existing regulations for laboratory safety regulations and good laboratory practice:



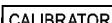
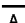
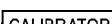
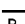

















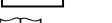
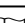
- First aid measures: In case of skin contact, immediately wash thoroughly with water and soap. Remove contaminated clothing and shoes and wash before reuse. If system fluid comes into contact with skin, wash thoroughly with water. After contact with the eyes carefully rinse the opened eye with running water for at least 10 minutes. Get medical attention if necessary.
- Personal precautions, protective equipment and emergency procedures:

Observe laboratory safety regulations. Avoid contact with skin and eyes. Do not swallow. Do not pipette by mouth. Do not eat, drink, smoke or apply makeup in areas where specimens or kit reagents are handled. When spilled, absorb with an inert material and put the spilled material in an appropriate waste disposal.

- Exposure controls / personal protection: Wear protective gloves of nitril rubber or natural latex. Wear protective glasses. Used according to intended use no dangerous reactions known.
- Conditions to avoid: Since substrate solution is light-sensitive. Store in the dark.
- For disposal of laboratory waste the national or regional legislation has to be observed.

Observe the guidelines for performing quality control in medical laboratories by assaying control sera.

## CONTENTS OF THE KIT

RUO 523	 96	Sufficient for 96 determinations
 MICROPLATE	1	One divisible microplate consisting of 12 modules of 8 wells each. Ready to use. Product code on module: <b>BPI</b>
 CALIBRATOR  A	1.5 ml	Calibrator A 0 U/ml, containing serum/buffer matrix (PBS, BSA, detergent, NaN3 0.09%), yellow. Ready to use.
 CALIBRATOR  B	1.5 ml	Calibrator B 6.3 U/ml, containing BPI antibodies in a serum/buffer matrix (PBS, BSA, detergent, NaN3 0.09%), yellow. Ready to use.
 CALIBRATOR  C	1.5 ml	Calibrator C 12.5 U/ml, containing BPI antibodies in a serum/buffer matrix (PBS, BSA, detergent, NaN3 0.09%), yellow. Ready to use.
 CALIBRATOR  D	1.5 ml	Calibrator D 25 U/ml, containing BPI antibodies in a serum/buffer matrix (PBS, BSA, detergent, NaN3 0.09%), yellow. Ready to use.
 CALIBRATOR  E	1.5 ml	Calibrator E 50 U/ml, containing BPI antibodies in a serum/buffer matrix (PBS, BSA, NaN3 0.09%), yellow. Ready to use.
 CALIBRATOR  F	1.5 ml	Calibrator F 100 U/ml, containing BPI antibodies in a serum/buffer matrix (PBS, BSA, detergent, NaN3 0.09%), yellow. Ready to use.
 CONTROL  +	1.5 ml	Control positive, containing BPI antibodies in a serum/buffer matrix (PBS, BSA, detergent, NaN3 0.09%), yellow. Ready to use. The concentration is specified on the certificate of analysis.
 CONTROL  -	1.5 ml	Control negative, containing BPI antibodies in a serum/buffer matrix (PBS, BSA, detergent, NaN3 0.09%), yellow. Ready to use. The concentration is specified on the certificate of analysis.
 DILUENT	20 ml	Sample Buffer P, containing PBS, BSA, detergent, preservative sodium azide 0.09%, yellow, concentrate (5 x).
 CONJUGATE	15 ml	Enzyme Conjugate containing anti-human IgG antibodies, HRP labelled; PBS, BSA, detergent, preservative PROCLIN 0.05%, light red. Ready to use.
 TMB	15 ml	TMB Substrate; containing 3,3', 5,5'- Tetramethylbenzidin, colorless. Ready to use.
 STOP	15 ml	Stop solution; contains acid. Ready to use.
 WASH	20 ml	Wash Buffer, containing Tris, detergent, preservative sodium azide 0.09%; 50 x conc.
 	1	Certificate of Analysis

## MATERIALS REQUIRED

- Microplate reader capable of endpoint measurements at 450 nm; optional: reference filter at 620 nm
- Data reduction software
- Multi-channel dispenser or repeatable pipette for 100 µl
- Vortex mixer
- Pipettes for 10 µl, 100 µl and 1000 µl
- Laboratory timing device
- Distilled or deionised water
- Measuring cylinder for 1000 ml and 100 ml
- Plastic container for storage of the wash solution

This ELISA assay is suitable for use on open automated ELISA processors. Each assay has to be validated on the respective automated system. Detailed information is provided upon request.

## SPECIMEN COLLECTION, STORAGE AND HANDLING

- Collect whole blood specimens using acceptable medical techniques to avoid hemolysis.
- Allow blood to clot and separate the serum or plasma by centrifugation.
- Test serum should be clear and non-hemolyzed. Contamination by hemolysis or lipemia should be avoided, but does not interfere with this assay.
- Specimens may be refrigerated at 2-8°C for up to five days or stored at -20°C up to six months.
- Avoid repetitive freezing and thawing of serum or plasma samples. This may result in variable loss of antibody activity.
- Testing of heat-inactivated sera is not recommended.

## STORAGE AND STABILITY

- Store test kit at 2-8°C in the dark.
- Do not expose reagents to heat, sun, or strong light during storage and usage.
- Store microplate sealed and dessicated in the clip bag provided.
- Shelf life of the unopened test kit is 18 months from day of production.  
Unopened reagents are stable until expiration of the kit. See labels for individual batch.
- Diluted Wash Buffer and Sample Buffer are stable for at least 30 days when stored at 2-8°C.  
We recommend consumption on the same day.

## PROCEDURAL NOTES

- Do not use kit components beyond their expiration dates.
- Do not interchange kit components from different lots and products.
- All materials must be at room temperature (20-28°C) prior to use.
- Prepare all reagents and samples. Once started, perform the test without interruption.
- Double determinations may be done. By this means pipetting errors may become obvious.
- Perform the assay steps only in the order indicated.
- Always use fresh sample dilutions.
- Pipette all reagents and samples into the bottom of the wells.
- To avoid carryover or contamination, change the pipette tip between samples and different kit controls.
- Wash microwells thoroughly and remove the last droplets of wash buffer.
- All incubation steps must be accurately timed.
- Do not re-use microplate wells.

## PREPARATION OF REAGENTS

### WASH

Dilute the contents of one vial of the buffered wash solution concentrate (50x) with distilled or deionised water to a final volume of 1000 ml prior to use.

### DILUENT

Sample Buffer P: Prior to use dilute the contents (20 ml) of one vial of sample buffer 5x concentrate with distilled or deionised water to a final volume of 100 ml.

### Preparation of samples

Dilute samples 1:100 before the assay: Put 990 µl of prediluted sample buffer in a polystyrene tube and add 10 µl of sample. Mix well. Note: Calibrators / Controls are ready to use and need not be diluted.

## TEST PROCEDURE

Prepare enough microplate modules for all calibrators / controls and samples.

1. Pipette 100 µl of calibrators, controls and prediluted samples into the wells.  
Incubate for 30 minutes at room temperature (20-28 °C).  
Discard the contents of the microwells and wash 3 times with 300 µl of wash solution.
2. Dispense 100 µl of enzyme conjugate into each well.  
Incubate for 15 minutes at room temperature.  
Discard the contents of the microwells and wash 3 times with 300 µl of wash solution.
3. Dispense 100 µl of TMB substrate solution into each well.  
Incubate for 15 minutes at room temperature
4. Add 100 µl of stop solution to each well of the modules  
Incubate for 5 minutes at room temperature.  
Read the optical density at 450 nm (reference 600-690nm) and calculate the results.  
The developed colour is stable for at least 30 minutes. Read during this time.

Example for a pipetting scheme:

	1	2	3	4	5	6	7	8	9	10	11	12
A	A	P1										
B	B	P2										
C	C	P3										
D	D											
E	E											
F	F											
G	C+											
H	C-											

P1, ... sample A-F calibrators C+, C- controls

## VALIDATION

Test results are valid if the optical densities at 450 nm for calibrators / controls and the results for controls comply with the reference ranges indicated on the Certificate of Analysis enclosed in each test kit.  
If these quality control criteria are not met the assay run is invalid and should be repeated.

## CALCULATION OF RESULTS

For quantitative results plot the optical density of each calibrator versus the calibrator concentration to create a calibration curve. The concentration of samples may then be estimated from the calibration curve by interpolation.  
Using data reduction software a 4-Parameter-Fit with lin-log coordinates for optical density and concentration is the data reduction method of choice.

## ASSAY CHARACTERISTICS

### Calibration

This assay system is calibrated in relative arbitrary units, since no international reference preparation is available for this assay.

### Measuring range

The calculation range of this ELISA assay is 0 - 100 U/ml

### Interfering substances

No interference has been observed with haemolytic (up to 1000 mg/dl) or lipemic (up to 3 g/dl triglycerides) sera or plasma, or bilirubin (up to 40 mg/dl) containing sera or plasma. Nor have any interfering effects been observed with the use of anticoagulants (Citrate, EDTA, Heparine). However for practical reasons it is recommended that grossly hemolyzed or lipemic samples should be avoided.

## LIMITATIONS OF THE PROCEDURE

For research use only. Not for diagnostic procedures.

## REFERENCES

1. Hauschild, S.; Schmitt, W.H.; Csernok, E.; Flesch, B. K.; Rautmann, A.; Gross, W.L. ANCA in systemic vasculitides, collagen vascular diseases, rheumatic disorders and inflammatory bowel diseases. ANCA-associated Vasculitides: Immunological and clinical aspects. Edited by W. L. Gross, Plenum Press, New York, 1993, pp. 245 - 251.
2. Schmitt, W.H.; Csernok, E.; Flesch, B. K.; Hauschild, S.; Gross, W.L. Autoantibodies directed against Lysozyme: A new target Antigen for anti-neutrophil cytoplasmic antibodies (ANCA). ANCA-associated Vasculitides: Immunological and clinical aspects. Edited by W. L. Gross, Plenum Press, New York, 1993, 267 - 272.
3. Skogh, T. and Peen, E. Lactoferrin, Anti-Lactoferrin antibodies and inflammatory disease. ANCA-associated Vasculitides: Immunological and clinical aspects. Edited by W. L. Gross, Plenum Press, New York, 1993, pp. 533 - 538.
4. Stoffel, M. P.; Csernok, E.; Herzberg, C.; Johnston, T.; Carroll, S.F. and Gross, W. L. Anti-neutrophil cytoplasmic antibodies (ANCA) directed against bactericidal/permeability increasing protein (BPI): a new seromarker for inflammatory bowel disease and associated disorders. Clin. Exp. Immunol. 1996, No. 104, pp. 54 - 59.
5. Zhao, M. H.; Jayne, D. R. W.; Ardiles, L. G.; Culley, F.; Hodson, M.E. and Lockwood. Autoantibodies against bactericidal/permeability increasing protein in patients with cystic fibrosis. J. Med. 1996, No. 89, pp. 259 - 265.
6. Wiik, A.; Stumman, L.; Kjeldsen, L.; Borregaard, N.; Ullman, S.; Jacobsen, S. and Halberg, P. The diversity of perinuclear antineutrophil cytoplasmic antibodies (pANCA) antigens. Proceedings of the 6th international ANCA workshop, France - July 1995, pp. 15 - 17.
7. Jenette, J. C. and Falk, R. J. Antineutrophil cytoplasmic autoantibodies and associated disease: A review. Am. J. of Kidney Disease 1990, Vol. XV, No. 6, pp. 517 - 529.
8. Schmitt, W. H. and Gross, W.L. Antineutrophil zytoplasmatische Antikörper (ANCA). Internist 1995, No. 36, pp. 282 - 290.
9. Gross, W.L. and Rheinhold-Keller, E. ANCA-assoziierte Vaskulitiden (Wegener Granulomatose, Churg-Strauss-Syndrom, Mikroskopische Polyangiitis) – 1. Systematik, Pathogenese und Klinik. Z. Rheumatol (1995), No. 54, pp. 279 - 290.
10. de Grot, K.; Schnabel, A. and Gross, W. L. ANCA-assoziierte Vaskulitiden (Wegener Granulomatose, Churg-Strauss-Syndrom, Mikroskopische Polyangiitis) - 1. Diagnostisches Procedere. Z. Rheumatol (1995), No. 54, pp. 291 - 302.
11. Schultz H, Weiss JP: The bactericidal/permeability-increasing protein (BPI) in infection and inflammatory disease. Clin Chim Acta 2007, 384:12-23.
12. Schultz H: From infection to autoimmunity: a new model for induction of ANCA against the bactericidal/permeability increasing protein (BPI). Autoimmun.Rev 2007, 6:223-227.
13. Talor MV, Stone JH, Stebbing J, Barin J, Rose NR, Burek CL: Antibodies to selected minor target antigens in patients with anti-neutrophil cytoplasmic antibodies (ANCA). Clin.Exp.Immunol. 2007, 150:42-48.

Change Control

Former version: -

Reason for revision: *first edition*