




Specificity Test (Cross Reaction Test with β -Klotho)

Report Made: 10th Dec 2020
 Code No. 27998
 Product Name: Human soluble α -Klotho Assay Kit - IBL

α -Klotho		β -Klotho	
pg/mL	A450	pg/mL	A450
6,000	2.391	6,000,000	0.006
3,000	1.444	3,000,000	-0.001
1,500	0.749	1,500,000	0.000
750	0.380	750,000	-0.003
375	0.185	375,000	-0.005
188	0.088	187,500	-0.002
94	0.045	93,750	-0.003

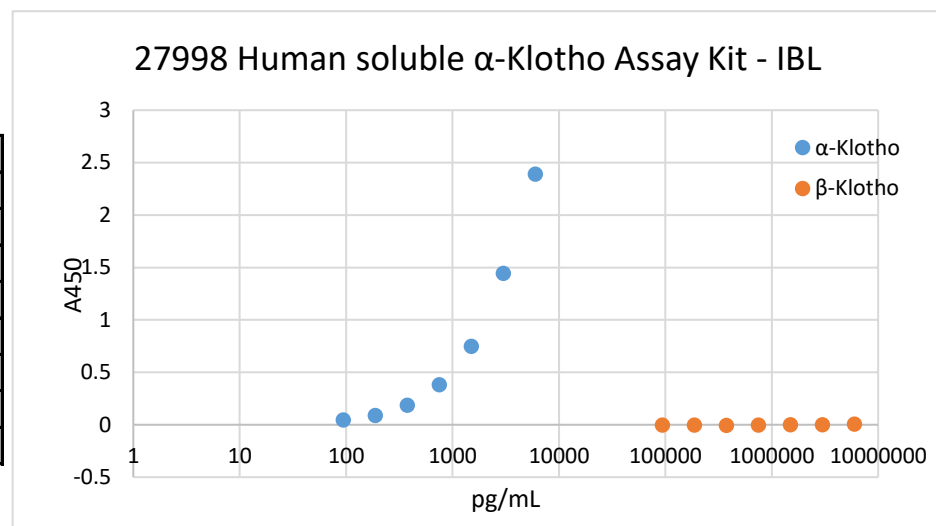

Distributed By:
IBL-America, Inc.
 8201 Central Ave NE, Suite P
 Minneapolis, MN 55432, USA
info@ibl-america.com
 (888) 523 1246