


Food Load test

	No.1	No.2	No.3	No.4	No.5	No.6
0hr	16.1	13.6	28.1	13.4	10.4	12.9
0.5hr	52.9	46.1	145.9	100.3	50.6	95.5
1hr	27.2	38.7	107.2	74.8	23.0	71.7
3hr	8.2	27.0	54.4	28.7	9.5	16.6

(pmol/L)

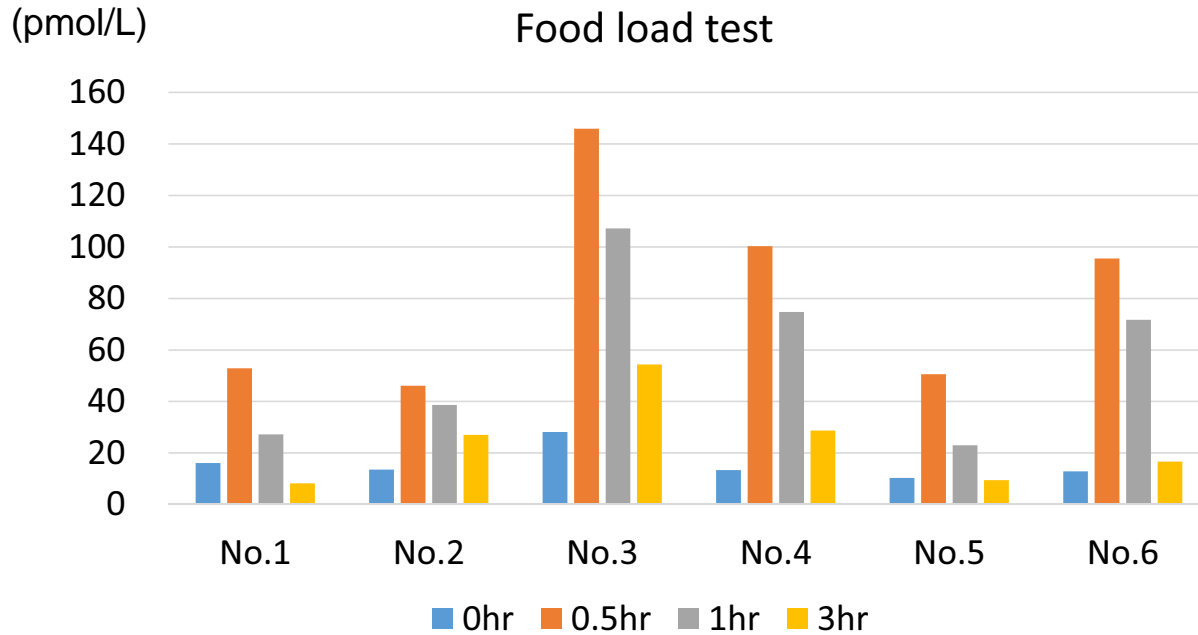

27702: Mouse GIP, Active
(High Sensitive) Assay Kit - IBL

<specimen>

ICR Mouse EDTA-Plasma
DPP-4 inhibitor(+)/Aprotinin(+)
(sampling from caudal vein)

<protocol>

- +
- Fasting for 16hr
- +
- Sampling (0hr)
- +
- Feeding
- +
- Sampling (0.5hr) and withdraw food
- +
- Sampling (1hr)
- +
- Sampling (3hr)

<sample volume>

10µL for duplicate sample examination
(20-fold dilution)

Distributed by:



Immuno-Biological Laboratories, Inc. (IBL-America)

8201 Central Ave NE, Suite P, Minneapolis, MN 55432, USA

Call Us : 888-523-1246 Email: info@ibl-america.com

URL: www.ibl-america.com