



# Human Bcl-2 ELISA

Enzyme-linked Immunosorbent Assay for  
the determination of human Bcl-2.

For research use only, not for use in diagnostic procedures.

**REF** *IB49641*

---

Manufactured for:  
**Immuno-Biological Laboratories, Inc. (IBL-America)**  
8201 Central Ave NE, Suite P, Minneapolis, MN 55432  
Toll Free: (888) 523-1246 Fax: (763) 780-2988  
[www.ibl-america.com](http://www.ibl-america.com) / [info@ibl-america.com](mailto:info@ibl-america.com)

## TABLE OF CONTENTS

1	Intended Use	3
2	Background	3
3	Principles of the Test	5
4	Reagents Provided	7
5	Storage Instructions – ELISA Kit	8
6	Specimen Collection and Storage Instructions	8
7	Materials Required But Not Provided	10
8	Precautions for Use	11
9	Preparation of Reagents	13
10	Test Protocol	17
11	Calculation of Results	22
12	Limitations	25
13	Performance Characteristics	26
14	Ordering Information	31
15	Reagent Preparation Summary	32
16	Test Protocol Summary	33

## 1 Intended Use

The human Bcl-2 ELISA is an enzyme-linked immunosorbent assay for the determination of human Bcl-2. **For research use only, not for use in diagnostic procedures.**

## 2 Background

Bcl-2 initiates a new gene family involved in the regulation of cell death and survival without affecting cell proliferation. This precise regulation and maintenance of balance between cell proliferation and cell death in multicellular organisms is critical for tissue homeostasis. Bcl-2, encoded by a proto-oncogene, is an intracellular membrane-associated protein that functions to block programmed cell death. Expression of Bcl-2 has been reported in a wide range of hematopoietic cells, nonneoplastic epithelia and epithelial malignancies. Bcl-2 levels are low or undetectable in endothelial cells.

The Bcl-2 gene has a unique function among mammalian oncogenes as a negative regulator of apoptosis. It was first discovered because of its involvement in chromosomal translocations commonly found in lymphomas . Bcl-2 is furthermore associated with stem cells committed to differentiation and morphogenesis.

A decrease in Bcl-2 levels leads to cell death by apoptosis. Overexpression of Bcl-2 on the other hand protects cells from death, but it is neither able to immortalize normal cells, nor to cause tumorigenic transformation of immortalized cells.

Heterogeneous expression of Bcl-2 in malignancies suggests that the gene is differentially regulated. Furthermore, its expression in association with precancerous lesions suggests a role in the early stage of tumorigenesis.

Several homologs of Bcl-2 have recently been discovered, some of which function as inhibitors of cell death and others as promoters of apoptosis that oppose the function of the Bcl-2 protein.

Bcl-2 protein plays a critical role in oncogenesis as well as resistance to therapy of cancer. Furthermore prognostic significance of Bcl-2

expression has been shown for several malignancies such as Non-Hodgkin's Lymphoma, squamous cell carcinomas, breast carcinoma, carcinoma of the gallbladder, and thymoma. Dysregulated Bcl-2 expression has been shown in Multiple Myeloma subjects and in acute Myeloid leukemia. Bcl-2 has been suggested as a useful marker for adequate IL-2 therapy e.g. in AIDS subjects.

For literature update refer to **[www.ibl-america.com](http://www.ibl-america.com)**

### 3 Principles of the Test

An anti-human Bcl-2 coating antibody is adsorbed onto microwells.

Human Bcl-2 present in the sample or standard binds to antibodies adsorbed to the microwells. A biotin-conjugated anti-human Bcl-2 antibody is added and binds to human Bcl-2 captured by the first antibody.

Following incubation unbound biotin-conjugated anti-human Bcl-2 antibody is removed during a wash step. Streptavidin-HRP is added and binds to the biotin-conjugated anti-human Bcl-2 antibody.

Following incubation unbound Streptavidin-HRP is removed during a wash step, and substrate solution reactive with HRP is added to the wells.

Figure 1

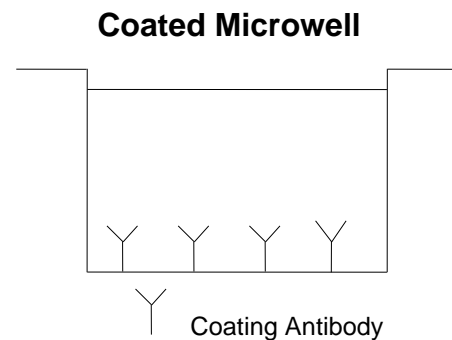


Figure 2

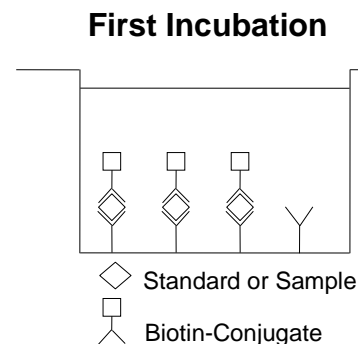


Figure 3



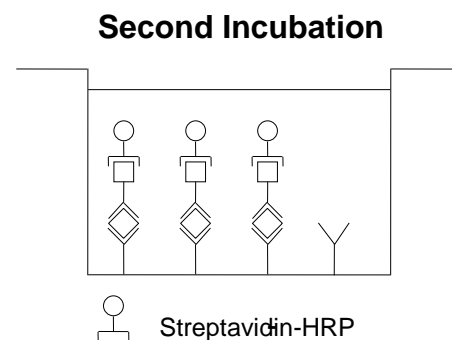


Figure 4

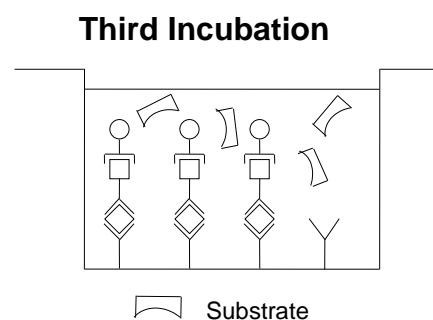

A coloured product is formed in proportion to the amount of human Bcl-2 present in the sample or standard. The reaction is terminated by addition of acid and absorbance is measured at 450 nm. A standard curve is prepared from 7 human Bcl-2 standard dilutions and human Bcl-2 sample concentration determined.

Figure 5

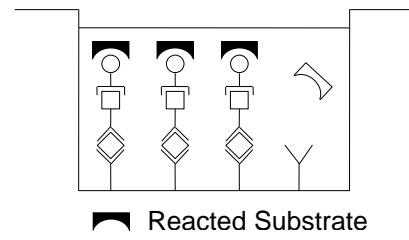

## 4 Reagents Provided

### 4.1 Reagents for human Bcl-2 ELISA IB49641 (96 tests)

- 1 aluminum pouch with a **Microwell Plate coated** with monoclonal antibody to human Bcl-2
- 1 vial (70 µl) **Biotin-Conjugate** anti-human Bcl-2 monoclonal antibody
- 1 vial (150 µl) **Streptavidin-HRP**
- 2 vials human Bcl-2 **Standard** lyophilized, 64 ng/ml upon reconstitution
- 1 vial (12 ml) **Sample Diluent**
- 1 vial (5 ml) **Assay Buffer Concentrate** 20x  
(PBS with 1% Tween 20, 10% BSA)
- 1 bottle (50 ml) **Wash Buffer Concentrate** 20x  
(PBS with 1% Tween 20)
- 1 bottle (15 ml) **Lysis Buffer** 10x
- 1 vial (15 ml) **Substrate Solution** (tetramethyl-benzidine)
- 1 vial (15 ml) **Stop Solution** (1M Phosphoric acid)
- 4 **Adhesive Films**

## 5 Storage Instructions – ELISA Kit

Store kit reagents between 2° and 8°C. Immediately after use remaining reagents should be returned to cold storage (2° to 8°C). Expiry of the kit and reagents is stated on labels.

Expiry of the kit components can only be guaranteed if the components are stored properly, and if, in case of repeated use of one component, this reagent is not contaminated by the first handling.

## 6 Specimen Collection and Storage Instructions

Cell lysates, serum and plasma were tested with this assay. Other biological samples might be suitable for use in the assay. Remove serum or plasma from the clot or cells as soon as possible after clotting and separation.

Samples containing a visible precipitate must be clarified prior to use in the assay. Do not use grossly hemolyzed or lipemic specimens.

Samples should be aliquoted and must be stored frozen at -20°C to avoid loss of bioactive human Bcl-2. If samples are to be run within 24 hours, they may be stored at 2° to 8°C (for sample stability refer to 13.5).

Avoid repeated freeze-thaw cycles. Prior to assay, the frozen sample should be brought to room temperature slowly and mixed gently.

### 6.1 Sample Preparation – Cell Lysate Protocol

Numerous extraction protocols can be used. The following protocol is provided as an example of a suitable extraction procedure, but should not be construed as necessarily being the method of choice. Users may wish to experiment with extraction procedures that work best in their hands.

For suspension cells: Pellet by centrifugation, remove supernatant and proceed to Addition of Lysis Buffer.

For attached cells: Remove supernatant from cells, wash cells once with PBS, harvest cells by scraping and gentle centrifugation, aspirate PBS, leaving an intact cell pellet (at this point the cell pellet can be frozen at -



80°C and lysed at a later date) and proceed to Addition of Lysis Buffer. Addition of Lysis Buffer: Resuspend the pellet in Lysis Buffer (1x) (we recommend a concentration of  $5 \times 10^6$  cells/ml.), incubate 60 minutes at room temperature with gentle shaking and transfer extracts to microcentrifuge tubes and centrifuge at 1000 x g for 15 minutes. Aliquot the cleared lysate to clean microfuge tubes and continue the test procedure (Alternatively lysates can be stored at -80°C and assayed at a later time. Divide lysates into small aliquots to avoid multiple freeze/thaw cycles.).

## **7 Materials Required But Not Provided**

- 5 ml and 10 ml graduated pipettes
- 5  $\mu$ l to 1000  $\mu$ l adjustable single channel micropipettes with disposable tips
- 50  $\mu$ l to 300  $\mu$ l adjustable multichannel micropipette with disposable tips
- Multichannel micropipette reservoir
- Beakers, flasks, cylinders necessary for preparation of reagents
- Device for delivery of wash solution (multichannel wash bottle or automatic wash system)
- Microwell strip reader capable of reading at 450 nm (620 nm as optional reference wave length)
- Glass-distilled or deionized water
- Statistical calculator with program to perform regression analysis

## 8 Precautions for Use

- All chemicals should be considered as potentially hazardous. We therefore recommend that this product is handled only by those persons who have been trained in laboratory techniques and that it is used in accordance with the principles of good laboratory practice. Wear suitable protective clothing such as laboratory overalls, safety glasses and gloves. Care should be taken to avoid contact with skin or eyes. In the case of contact with skin or eyes wash immediately with water. See material safety data sheet(s) and/or safety statement(s) for specific advice.
- Reagents are intended for research use only, not for use in diagnostic procedures.
- Do not mix or substitute reagents with those from other lots or other sources.
- Do not use kit reagents beyond expiration date on label.
- Do not expose kit reagents to strong light during storage or incubation.
- Do not pipette by mouth.
- Do not eat or smoke in areas where kit reagents or samples are handled.
- Avoid contact of skin or mucous membranes with kit reagents or specimens.
- Rubber or disposable latex gloves should be worn while handling kit reagents or specimens.
- Avoid contact of substrate solution with oxidizing agents and metal.
- Avoid splashing or generation of aerosols.
- In order to avoid microbial contamination or cross-contamination of reagents or specimens which may invalidate the test use disposable pipette tips and/or pipettes.
- Use clean, dedicated reagent trays for dispensing the conjugate and substrate reagent.

- Exposure to acid inactivates the conjugate.
- Glass-distilled water or deionized water must be used for reagent preparation.
- Substrate solution must be at room temperature prior to use.
- Decontaminate and dispose specimens and all potentially contaminated materials as they could contain infectious agents. The preferred method of decontamination is autoclaving for a minimum of 1 hour at 121.5°C.
- Liquid wastes not containing acid and neutralized waste may be mixed with sodium hypochlorite in volumes such that the final mixture contains 1.0% sodium hypochlorite. Allow 30 minutes for effective decontamination. Liquid waste containing acid must be neutralized prior to the addition of sodium hypochlorite.

## 9 Preparation of Reagents

**Buffer concentrates** should be brought to room temperature and should be diluted before starting the test procedure.

If crystals have formed in the **Buffer Concentrates**, warm them gently until they have completely dissolved.

### 9.1 Wash Buffer (1x)

Pour entire contents (50 ml) of the **Wash Buffer Concentrate (20x)** into a clean 1000 ml graduated cylinder. Bring to final volume of 1000 ml with glass-distilled or deionized water.

Mix gently to avoid foaming.

Transfer to a clean wash bottle and store at 2° to 25°C. Please note that Wash Buffer (1x) is stable for 30 days.

Wash Buffer (1x) may also be prepared as needed according to the following table:

Number of Strips	Wash Buffer Concentrate (20x) (ml)	Distilled Water (ml)
1 - 6	25	475
1 - 12	50	950

### 9.2 Assay Buffer(1x)

Pour the entire contents (5 ml) of the **Assay Buffer Concentrate(20x)** into a clean 100 ml graduated cylinder. Bring to final volume of 100 ml with distilled water. Mix gently to avoid foaming.

Store at 2° to 8°C. Please note that the Assay Buffer(1x) is stable for 30 days.

Assay Buffer(1x) may also be prepared as needed according to the following table:

Number of Strips	Assay Buffer Concentrate (20x) (ml)	Distilled Water (ml)
1 - 6	2.5	47.5
1 - 12	5.0	95.0

### 9.3 Lysis Buffer

Pour the entire contents (15 ml) of the **Lysis Buffer Concentrate** (10x) into a clean 150 ml graduated cylinder. Bring to final volume of 150 ml with distilled or deionized water and mix gently. Store at room temperature. Please note that the Lysis Buffer is stable for 30 days.

### 9.4 Biotin-Conjugate

**Please note that the Biotin-Conjugate should be used within 30 minutes after dilution.**

Make a 1:100 dilution of the concentrated **Biotin-Conjugate** solution with Assay Buffer (1x) in a clean plastic tube as needed according to the following table:

Number of Strips	Biotin-Conjugate (ml)	Assay Buffer Red (1x) (ml)
1 - 6	0.03	2.97
1 - 12	0.06	5.94

## 9.5 Streptavidin-HRP

**Please note that the Streptavidin-HRP should be used within 30 minutes after dilution.**

Make a 1:100 dilution of the concentrated **Streptavidin-HRP** solution with Assay Buffer Red (1x) in a clean plastic tube as needed according to the following table:

Number of Strips	Streptavidin-HRP (ml)	Assay Buffer Red (1x) (ml)
1 - 6	0.06	5.94
1 - 12	0.12	11.88

## 9.6 Human Bcl-2 Standard

Reconstitute **human Bcl-2 standard** by addition of distilled water. Reconstitution volume is stated on the label of the standard vial. Swirl or mix gently to insure complete and homogeneous solubilization (concentration of reconstituted standard = 64 ng/ml). Allow the standard to reconstitute for 10-30 minutes. Mix well prior to making dilutions.

After usage remaining standard cannot be stored and has to be discarded.

**Standard dilutions** can be prepared directly on the microwell plate (see 10.d) or alternatively in tubes (see 9.6.1).

### 9.6.1 External Standard Dilution

Label 7 tubes, one for each standard point.

S1, S2, S3, S4, S5, S6, S7

Then prepare 1:2 serial dilutions for the standard curve as follows:  
Pipette 225 µl of Sample Diluent into each tube.


Pipette 225 µl of reconstituted standard (concentration = 64 ng/ml) into the first tube, labelled S1, and mix (concentration of standard 1 = 32 ng/ml).

Pipette 225  $\mu$ l of this dilution into the second tube, labelled S2, and mix thoroughly before the next transfer.

Repeat serial dilutions 5 more times thus creating the points of the standard curve (see Figure 6).

Sample Diluent serves as blank.

Figure 6

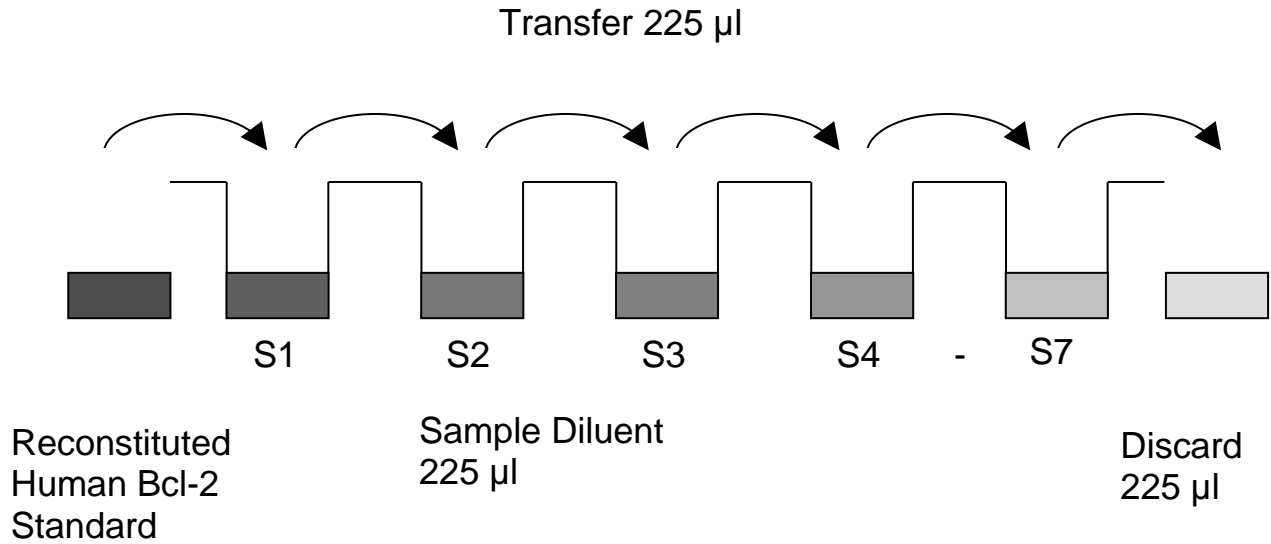




## 10 Test Protocol

- a. For cell lysis follow the cell lysate protocol (see 6.1 Sample Preparation – Cell Lysate Protocol):
- b. Determine the number of microwell strips required to test the desired number of samples plus appropriate number of wells needed for running blanks and standards. Each sample, standard, blank and optional control sample should be assayed in duplicate. Remove extra microwell strips from holder and store in foil bag with the desiccant provided at 2°-8°C sealed tightly.
- c. Wash the microwell strips twice with approximately 400 µl **Wash Buffer** per well with thorough aspiration of microwell contents between washes. Allow the Wash Buffer to sit in the wells for about **10 – 15 seconds** before aspiration. Take care not to scratch the surface of the microwells.  
After the last wash step, empty wells and tap microwell strips on absorbent pad or paper towel to remove excess Wash Buffer. Use the microwell strips immediately after washing. Alternatively microwell strips can be placed upside down on a wet absorbent paper for not longer than 15 minutes. **Do not allow wells to dry.**
- d. **Standard dilution on the microwell plate** (Alternatively the standard dilution can be prepared in tubes - see 9.6.1):  
Add 100 µl of Sample Diluent in duplicate to all **standard wells**.  
Pipette 100 µl of prepared **standard** (see Preparation of Standard 9.6, concentration = 64 ng/ml) in duplicate into well A1 and A2 (see Table 1). Mix the contents of wells A1 and A2 by repeated aspiration and ejection (concentration of standard 1, S1 = 32 ng/ml), and transfer 100 µl to wells B1 and B2, respectively (see Figure 7). Take care not to scratch the inner surface of the microwells. Continue this procedure 5 times, creating two rows of human Bcl-2 standard dilutions ranging from 32.0 to 0.5 ng/ml. Discard 100 µl of the contents from the last microwells (G1, G2) used.

Figure 7

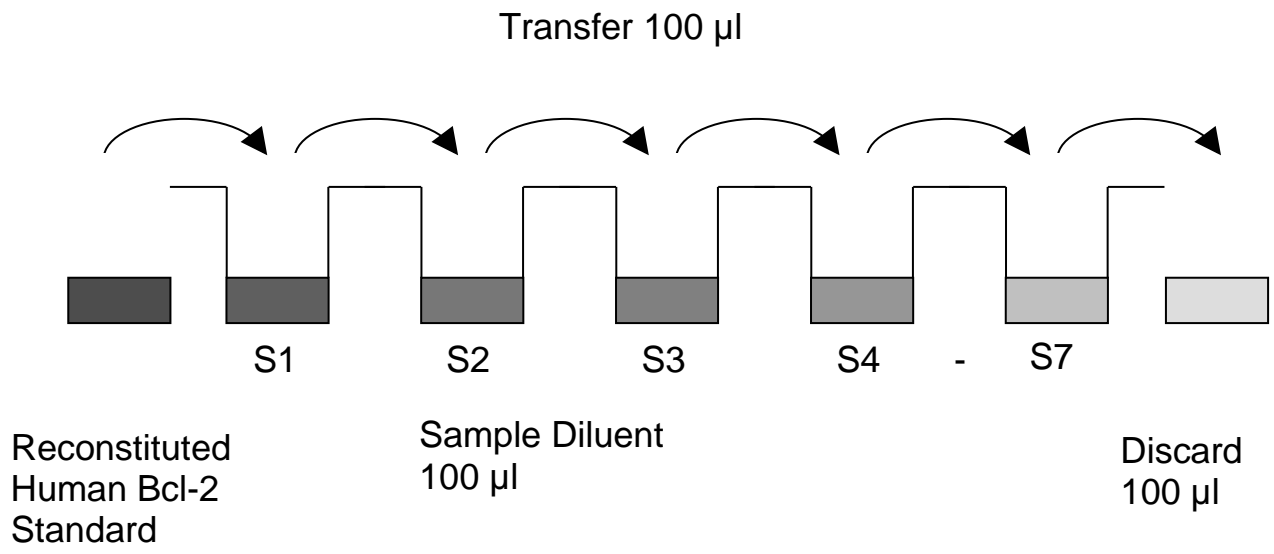

In case of an **external standard dilution** (see 9.6.1), pipette 100 µl of these standard dilutions (S1 - S7) in the standard wells according to Table 1.

Table 1

Table depicting an example of the arrangement of blanks, standards and samples in the microwell strips:

	<b>1</b>	<b>2</b>	<b>3</b>	<b>4</b>
<b>A</b>	Standard 1 (32.0 ng/ml)	Standard 1 (32.0 ng/ml)	Sample 1	Sample 1
<b>B</b>	Standard 2 (16.0 ng/ml)	Standard 2 (16.0 ng/ml)	Sample 2	Sample 2
<b>C</b>	Standard 3 (8.0 ng/ml)	Standard 3 (8.0 ng/ml)	Sample 3	Sample 3
<b>D</b>	Standard 4 (4.0 ng/ml)	Standard 4 (4.0 ng/ml)	Sample 4	Sample 4
<b>E</b>	Standard 5 (2.0 ng/ml)	Standard 5 (2.0 ng/ml)	Sample 5	Sample 5
<b>F</b>	Standard 6 (1.0 ng/ml)	Standard 6 (1.0 ng/ml)	Sample 6	Sample 6
<b>G</b>	Standard 7 (0.5 ng/ml)	Standard 7 (0.5 ng/ml)	Sample 7	Sample 7
<b>H</b>	Blank	Blank	Sample 8	Sample 8

- e. Add 100 µl of **Sample Diluent** in duplicate to the **blank wells**.
- f. Add 80 µl of **Sample Diluent** to the **sample wells**.
- g. Add 20 µl of each **sample** in duplicate to the **sample wells**.
- h. Prepare **Biotin-Conjugate** (see Preparation of Biotin-Conjugate 9.4).
- i. Add 50 µl of **Biotin-Conjugate** to all wells.
- j. Cover with an adhesive film and incubate at room temperature (18 to 25°C) for 2 hours, if available on a microplate shaker set at 400 rpm.
- k. Prepare **Streptavidin-HRP** (refer to Preparation of Streptavidin-HRP 9.5).
- l. Remove adhesive film and empty wells. **Wash** microwell strips 3 times according to point c. of the test protocol. Proceed immediately to the next step.
- m. Add 100 µl of diluted **Streptavidin-HRP** to all wells, including the blank wells.
- n. Cover with an adhesive film and incubate at room temperature (18° to 25°C) for 1 hour, if available on a microplate shaker set at 400 rpm.
- o. Remove adhesive film and empty wells. **Wash** microwell strips 3 times according to point c. of the test protocol. Proceed immediately to the next step.
- p. Pipette 100 µl of **TMB Substrate Solution** to all wells.
- q. Incubate the microwell strips at room temperature (18° to 25°C) for about 10 min. Avoid direct exposure to intense light.

**The colour development on the plate should be monitored and the substrate reaction stopped (see next point of this protocol) before positive wells are no longer properly recordable. Determination of the ideal time period for colour development has to be done individually for each assay.**

It is recommended to add the stop solution when the highest standard has developed a dark blue colour. Alternatively the colour development can be monitored by the ELISA reader at 620 nm. The substrate reaction should be stopped as soon as Standard 1 has reached an OD of 0.9 – 0.95.

- r. Stop the enzyme reaction by quickly pipetting 100 µl of **Stop Solution** into each well. It is important that the Stop Solution is spread quickly and uniformly throughout the microwells to completely inactivate the enzyme. Results must be read immediately after the Stop Solution is added or within one hour if the microwell strips are stored at 2 - 8°C in the dark.
- s. Read absorbance of each microwell on a spectro-photometer using 450 nm as the primary wave length (optionally 620 nm as the reference wave length; 610 nm to 650 nm is acceptable). Blank the plate reader according to the manufacturer's instructions by using the blank wells. Determine the absorbance of both the samples and the standards.

**Note: In case of incubation without shaking the obtained O.D. values may be lower than indicated below. Nevertheless the results are still valid.**

## 11 Calculation of Results

- Calculate the average absorbance values for each set of duplicate standards and samples. Duplicates should be within 20 per cent of the mean value.
- Create a standard curve by plotting the mean absorbance for each standard concentration on the ordinate against the human Bcl-2 concentration on the abscissa. Draw a best fit curve through the points of the graph (a 5-parameter curve fit is recommended).
- To determine the concentration of circulating human Bcl-2 for each sample, first find the mean absorbance value on the ordinate and extend a horizontal line to the standard curve. At the point of intersection, extend a vertical line to the abscissa and read the corresponding human Bcl-2 concentration.
- **If instructions in this protocol have been followed samples have been diluted 1:5 (20 µl sample + 80 µl Sample Diluent), the concentration read from the standard curve must be multiplied by the dilution factor (x 5).**
- **Calculation of samples with a concentration exceeding standard 1 may result in incorrect, low human Bcl-2 levels. Such samples require further external predilution according to expected human Bcl-2 values with Sample Diluent in order to precisely determine the actual human Bcl-2 level.**
- It is suggested that each testing facility establishes a control sample of known human Bcl-2 concentration and runs this additional control with each assay. If the values obtained are not within the expected range of the control, the assay results may be invalid.
- A representative standard curve is shown in Figure 8. This curve cannot be used to derive test results. Each laboratory must prepare a standard curve for each group of microwell strips assayed.

Figure 8

Representative standard curve for human Bcl-2 ELISA. Human Bcl-2 was diluted in serial 2-fold steps in Sample Diluent. Do not use this standard curve to derive test results. A standard curve must be run for each group of microwell strips assayed.

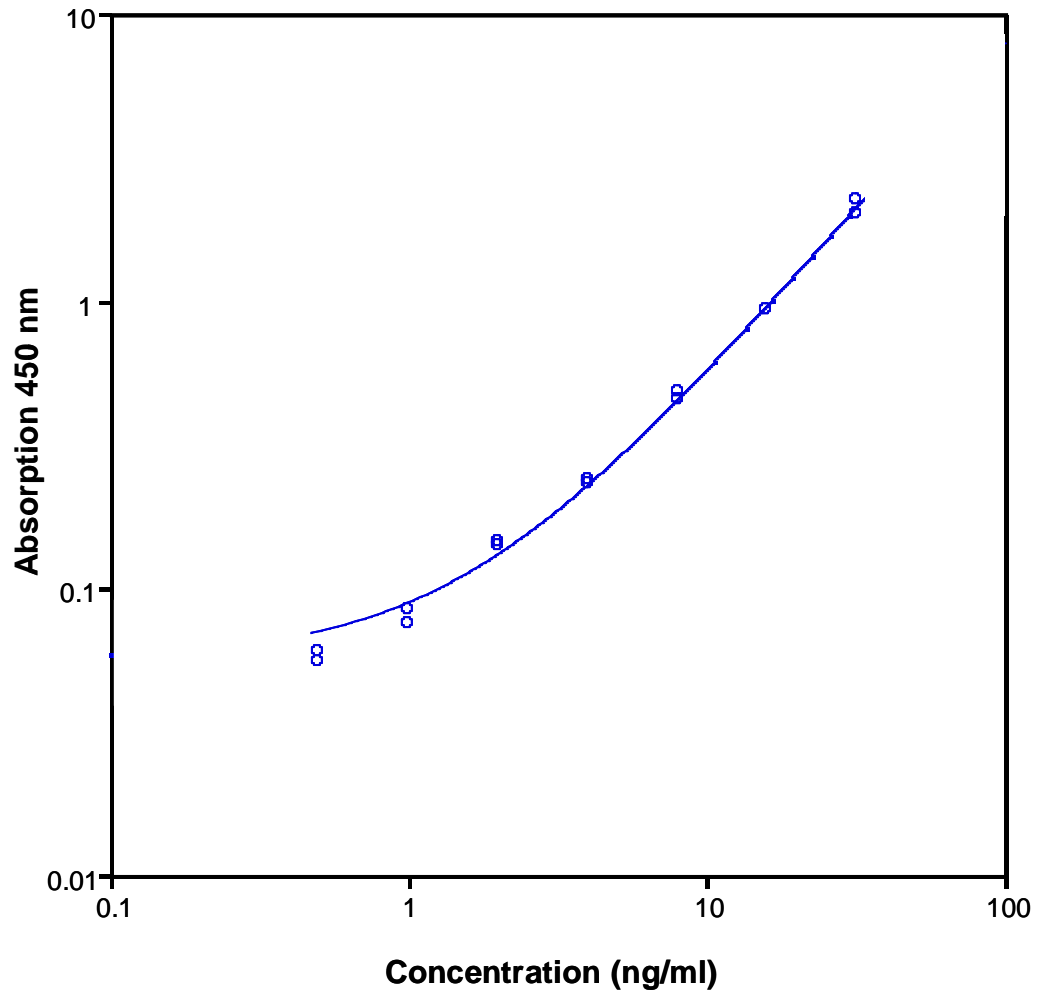


Table 2

Typical data using the human Bcl-2 ELISA

Measuring wavelength: 450 nm

Reference wavelength: 620 nm

Standard	Human Bcl-2 Concentration (ng/ml)	O.D. at 450 nm	Mean O.D. at 450 nm	C.V. (%)
1	32.0	2.248 2.007	2.128	5.7
2	16.0	0.937 0.940	0.939	0.2
3	8.0	0.453 0.483	0.468	3.2
4	4.0	0.231 0.240	0.236	1.9
5	2.0	0.141 0.145	0.143	1.4
6	1.0	0.075 0.084	0.080	5.7
7	0.5	0.060 0.056	0.058	3.4
Blank	0.0	0.018 0.018	0.018	0.0

The OD values of the standard curve may vary according to the conditions of assay performance (e.g. operator, pipetting technique, washing technique or temperature effects). Furthermore shelf life of the kit may affect enzymatic activity and thus colour intensity. Values measured are still valid.



## 12 Limitations

- Since exact conditions may vary from assay to assay, a standard curve must be established for every run.
- Bacterial or fungal contamination of either screen samples or reagents or cross-contamination between reagents may cause erroneous results.
- Disposable pipette tips, flasks or glassware are preferred, reusable glassware must be washed and thoroughly rinsed of all detergents before use.
- Improper or insufficient washing at any stage of the procedure will result in either false positive or false negative results. Empty wells completely before dispensing fresh wash solution, fill with Wash Buffer as indicated for each wash cycle and do not allow wells to sit uncovered or dry for extended periods.
- The use of radioimmunotherapy has significantly increased the number of subjects with human anti-mouse IgG antibodies (HAMA). HAMA may interfere with assays utilizing murine monoclonal antibodies leading to both false positive and false negative results. Serum samples containing antibodies to murine immunoglobulins can still be analysed in such assays when murine immunoglobulins (serum, ascitic fluid, or monoclonal antibodies of irrelevant specificity) are added to the sample.

## **13 Performance Characteristics**

### **13.1 Sensitivity**

The limit of detection of human Bcl-2 defined as the analyte concentration resulting in an absorbance significantly higher than that of the dilution medium (mean plus 2 standard deviations) was determined to be < 0.5 ng/ml (mean of 6 independent assays).

### **13.2 Reproducibility**

#### **13.2.1 Intra-assay**

Reproducibility within the assay was evaluated in 3 independent experiments. Each assay was carried out with 6 replicates of 8 samples containing different concentrations of human Bcl-2. 2 standard curves were run on each plate. Data below show the mean human Bcl-2 concentration and the coefficient of variation for each sample (see Table 3). The calculated overall intra-assay coefficient of variation was 8.6%.

Table 3

The mean human Bcl-2 concentration and the coefficient of variation for each sample

Sample	Experiment	Mean Human Bcl-2 Concentration (ng/ml)	Coefficient of Variation (%)
1	1	87.9	5.0
	2	76.6	9.0
	3	64.3	4.0
2	1	63.4	7.0
	2	54.4	9.0
	3	47.9	9.0
3	1	54.0	14.0
	2	40.8	9.0
	3	40.0	14.0
4	1	44.0	11.0
	2	35.3	7.0
	3	34.8	12.0
5	1	42.0	4.0
	2	46.1	4.0
	3	53.8	9.0
6	1	16.0	10.0
	2	19.8	7.0
	3	16.7	10.0
7	1	7.4	6.0
	2	7.7	11.0
	3	8.3	10.0
8	1	5.0	11.0
	2	5.5	7.0
	3	5.3	7.0

### 13.2.2 Inter-assay

Assay to assay reproducibility within one laboratory was evaluated in 3 independent experiments. Each assay was carried out with 6 replicates of 8 samples containing different concentrations of human Bcl-2. 2 standard curves were run on each plate. Data below show the mean human Bcl-2 concentration and the coefficient of variation calculated on 18 determinations of each sample (see Table 4). The calculated overall inter-assay coefficient of variation was 12.0%.

Table 4

The mean human Bcl-2 concentration and the coefficient of variation of each sample

Sample	Mean Human Bcl-2 Concentration (ng/ml)	Coefficient of Variation (%)
1	75.3	15.8
2	55.2	14.0
3	44.8	17.7
4	38.0	13.5
5	47.3	12.6
6	17.5	11.6
7	7.8	5.9
8	5.3	5.1

### 13.3 Spike Recovery

The spike recovery was evaluated by spiking 4 levels of human Bcl-2 into lysis buffer. Recoveries were determined in 3 independent experiments with 4 replicates each.

The overall mean recovery was 71%.

### 13.4 Dilution Parallelism

4 samples with different levels of human Bcl-2 were analysed at serial 2 fold dilutions with 4 replicates each.

The recovery ranged from 85% to 123% with an overall recovery of 102% (see Table 5).

Table 5

Sample	Dilution	Expected Human Bcl-2 Concentration (ng/ml)	Observed Human Bcl-2 Concentration (ng/ml)	Recovery of Expected Human Bcl-2 Concentration (%)
1	1:5	-	81.7	-
	1:10	40.9	34.5	85
	1:20	20.4	22.4	110
	1:40	10.2	11.3	111
2	1:5	-	50.3	-
	1:10	25.2	23.3	92
	1:20	12.6	12.8	102
	1:40	6.3	6.4	102
3	1:5	-	138.0	-
	1:10	69.0	75.6	110
	1:20	34.5	42.4	123
	1:40	17.3	19.6	113
4	1:5	-	113.5	-
	1:10	56.8	53.3	94
	1:20	28.4	25.7	91
	1:40	14.2	13.0	92

## **13.5 Sample Stability**

### **13.5.1 Freeze-Thaw Stability**

Aliquots of samples were stored at  $-20^{\circ}\text{C}$  and thawed 5 times, and the human Bcl-2 levels determined. There was no significant loss of human Bcl-2 immunoreactivity detected by freezing and thawing.

### **13.5.2 Storage Stability**

Aliquots of samples were stored at  $-20^{\circ}\text{C}$ ,  $2-8^{\circ}\text{C}$ , room temperature (RT) and at  $37^{\circ}\text{C}$ , and the human Bcl-2 level determined after 24 h. There was no significant loss of human Bcl-2 immunoreactivity detected during storage at  $-20^{\circ}\text{C}$ . A significant loss of human Bcl-2 immunoreactivity was detected during storage at  $4^{\circ}\text{C}$  (30%), at RT (45%) and  $37^{\circ}\text{C}$  (70%) after 24 h.

## **14 Ordering Information**

This kit is manufactured for Immuno-Biological Laboratories, Inc. (IBL-America). For ordering information, please contact:

**Immuno-Biological Laboratories, Inc. (IBL-America)**

8201 Central Ave NE, Suite P

Minneapolis, MN 55432

Toll Free: (888) 523-1246

Fax: (763) 780-2988

[www.ibl-america.com](http://www.ibl-america.com)

[info@ibl-america.com](mailto:info@ibl-america.com)

## 15 Reagent Preparation Summary

### 15.1 Wash Buffer (1x)

Add **Wash Buffer Concentrate 20x** (50 ml) to 950 ml distilled water.

Number of Strips	Wash Buffer Concentrate (ml)	Distilled Water (ml)
1 - 6	25	475
1 - 12	50	950

### 15.2 Assay Buffer(1x)

Add **Assay Buffer Concentrate 20x** (5 ml) to 95 ml distilled water.

Number of Strips	Assay Buffer Concentrate (ml)	Distilled Water (ml)
1 - 6	2.5	47.5
1 - 12	5.0	95.0

### 15.3 Lysis Buffer

Add **Lysis Buffer Concentrate 10x** (15 ml) to 135 ml distilled water.

### 15.4 Biotin-Conjugate

Make a 1:100 dilution of **Biotin-Conjugate** in Assay Buffer Red (1x):

Number of Strips	Biotin-Conjugate (ml)	Assay Buffer Red (1x) (ml)
1 - 6	0.03	2.97
1 - 12	0.06	5.94

### 15.5 Streptavidin-HRP

Make a 1:100 dilution of **Streptavidin-HRP** in Assay Buffer Red (1x):

Number of Strips	Streptavidin-HRP (ml)	Assay Buffer Red (1x) (ml)
1 - 6	0.06	5.94
1 - 12	0.12	11.88

### 15.6 Human Bcl-2 Standard

Reconstitute lyophilized **human Bcl-2 standard** with distilled water.  
(Reconstitution volume is stated on the label of the standard vial.)



## 16 Test Protocol Summary

1. Prepare cell extracts according to the cell lysis protocol.
2. Determine the number of microwell strips required.
3. Wash microwell strips twice with Wash Buffer.
4. Standard dilution on the microwell plate: Add 100 µl Sample Diluent, in duplicate, to all standard wells. Pipette 100 µl prepared standard into the first wells and create standard dilutions by transferring 100 µl from well to well. Discard 100 µl from the last wells.  
Alternatively external standard dilution in tubes (see 9.6.1): Pipette 100 µl of these standard dilutions in the microwell strips.
5. Add 100 µl Sample Diluent, in duplicate, to the blank wells.
6. Add 80 µl Sample Diluent to sample wells.
7. Add 20 µl sample in duplicate, to designated sample wells.
8. Prepare Biotin-Conjugate.
9. Add 50 µl Biotin-Conjugate to all wells.
10. Cover microwell strips and incubate 2 hours at room temperature (18° to 25°C).
11. Prepare Streptavidin-HRP.
12. Empty and wash microwell strips 3 times with Wash Buffer.
13. Add 100 µl diluted Streptavidin-HRP to all wells.
14. Cover microwell strips and incubate 1 hour at room temperature (18° to 25°C).
15. Empty and wash microwell strips 3 times with Wash Buffer.
16. Add 100 µl of TMB Substrate Solution to all wells.
17. Incubate the microwell strips for about 10 minutes at room temperature (18° to 25°C).
18. Add 100 µl Stop Solution to all wells.
19. Blank microwell reader and measure colour intensity at 450 nm.

**Note: If instructions in this protocol have been followed samples have been diluted 1:5 (20 µl sample + 80 µl Sample Diluent), the concentration read from the standard curve must be multiplied by the dilution factor (x 5).**