



## User's Manual

# Borrelia burgdorferi IgG ELISA

Enzyme immunoassay for the detection of human IgG antibodies against *Borrelia burgdorferi* in serum and plasma.



**IB79804**



**96 tests**

**RUO**

**For Research Use Only – Not for Use in Diagnostic Procedures**

## Contents

1. INTENDED USE .....	3
2. GENERAL INFORMATION .....	3
3. PRINCIPLE OF THE TEST .....	3
4. LIMITATIONS, PRECAUTIONS AND GENERAL COMMENTS .....	4
5. REAGENTS PROVIDED .....	4
6. MATERIALS REQUIRED BUT NOT PROVIDED .....	5
7. SAMPLE COLLECTION AND HANDLING .....	6
8. ASSAY PROCEDURE .....	6
9. EVALUATION .....	7
10. ASSAY CHARACTERISTICS .....	7
11. REFERENCES .....	8
SYMBOLS USED WITH IBL-AMERICA ASSAYS .....	9

## 1. INTENDED USE

The IBL-AMERICA *Borrelia burgdorferi* IgG antibody ELISA kit has been designed for the detection of specific IgG antibodies against *Borrelia burgdorferi* in serum and plasma. For research use only, not for use in diagnostic procedures.

## 2. GENERAL INFORMATION

*Borrelia burgdorferi* belongs to the family of spirochetes, of which three types have been identified to be human pathogenic: *Borrelia burgdorferi sensu stricto*, *B. afzelii* and *B. garinii*.

The *Borrelia burgdorferi* IgG ELISA test kit contains besides a whole cell antigen extract of *Borrelia burgdorferi sensu stricto*, which cross-reacts with *Borrelia afzelii* and *Borrelia garinii*, an addition of pure OspC, which increases the specificity and sensitivity of the assay.

## 3. PRINCIPLE OF THE TEST

The *Borrelia burgdorferi* IgG antibody test kit is based on the principle of the enzyme immunoassay (EIA). *Borrelia burgdorferi* antigen is bound on the surface of the microtiter strips. Diluted serum or ready-to-use standards are pipetted into the wells of the microtiter plate. A binding between the IgG antibodies of the serum and the immobilized *Borrelia burgdorferi* antigen takes place. After a one hour incubation at room temperature, the plate is rinsed with diluted wash solution, in order to remove unbound material. Then ready-to-use anti-human-IgG peroxidase conjugate is added and incubated for 30 minutes. After a further washing step, the substrate (TMB) solution is pipetted and incubated for 20 minutes, inducing the development of a blue dye in the wells. The color development is terminated by the addition of a stop solution, which changes the color from blue to yellow. The resulting dye is measured spectrophotometrically at the wavelength of 450 nm. The concentration of IgG antibodies is directly proportional to the intensity of the color.

#### 4. LIMITATIONS, PRECAUTIONS AND GENERAL COMMENTS

- This assay is intended for research use only – not for use in diagnostic procedures.
- Do not ingest or swallow! The usual laboratory safety precautions as well as the prohibition of eating, drinking and smoking in the lab have to be followed.
- All sera and plasma or buffers based upon, have been tested respective to HBsAg, HIV and HCV with recognized methods and were found negative. Nevertheless precautions like the use of latex gloves have to be taken.
- Serum and reagent spills have to be wiped off with a disinfecting solution (e.g. sodium hypochlorite, 5%) and have to be disposed of properly.
- All reagents have to be brought to room temperature (18 to 25 °C) before performing the test.
- Before pipetting all reagents should be mixed thoroughly by gentle tilting or swinging. Vigorous shaking with formation of foam should be avoided.
- It is important to pipet with constant intervals, so that all the wells of the microtiter plate have the same conditions.
- When removing reagents out of the bottles, care has to be taken that the stoppers are not contaminated. Further a possible mix-up has to be avoided. The content of the bottles is usually sensitive to oxidation, so that they should be opened only for a short time.
- In order to avoid a carry-over or a cross-contamination, separate disposable pipet tips have to be used.
- No reagents from different kit lots have to be used, they should not be mixed among one another.
- All reagents have to be used within the expiry period.
- In accordance with Good Laboratory Practice (GLP) or following ISO9001 all laboratory devices employed should be regularly checked regarding the accuracy and precision. This refers amongst others to microliter pipets and washing or reading (ELISA-Reader) instrumentation.
- The contact of certain reagents, above all the stopping solution and the substrate with skin, eye and mucosa has to be avoided, because possible irritations and acid burns could arise, and there exists a danger of intoxication.

#### 5. REAGENTS PROVIDED

Symbol	Components	Volume / Qty.
<b>SORB</b> <b>MT</b>	<i>Borrelia burgdorferi</i> antigen coated microtiter strips	12
<b>CAL</b> <b>A</b>	Calibrator A (Negative Control)	2 mL
<b>CAL</b> <b>B</b>	Calibrator B (Cut-Off Standard)	2 mL
<b>CAL</b> <b>C</b>	Calibrator C (Weak Positive Control)	2 mL
<b>CAL</b> <b>D</b>	Calibrator D (Positive Control)	2 mL
<b>ENZ</b> <b>CONJ</b>	Enzyme Conjugate	15 mL
<b>SUB</b> <b>TMB</b>	Substrate	15 mL
<b>STOP</b> <b>SOLN</b>	Stop Solution	15 mL
<b>SAM</b> <b>DIL</b>	Sample Diluent	60 mL
<b>WASH</b> <b>SOLN</b> <b>10x</b>	Washing Buffer (10x)	60 mL

#### Storage and Stability (refer to the expiry date on the outer box label)

Store kit components at 2-8°C and do not use after the expiry date on the box outer label. Before use, all components should be allowed to warm up to ambient temperature (18-25°C). After use, the plate should be resealed, the bottle caps replaced and tightened and the kit stored at 2-8°C. After the first opening the kit should be used within 3 months, the diluted wash buffer can be kept for 4 weeks at 2-8°C.

**5.1. Microtiter Strips**

12 strips with 8 breakable wells each, coated with *Borrelia* antigen (mix of *Borrelia sensu stricto*, *Borrelia afzelii* und *Borrelia garinii*). Ready-to-use

**5.2. Calibrator A (Negative Control)**

2 mL, protein solution diluted with PBS, contains no IgG antibodies against *Borrelia burgdorferi*. Addition of 0.01 % methylisothiazolone and 0.01 % bromonitrodioxane. Ready-to-use.

**5.3. Calibrator B (Cut-Off Standard)**

2 mL human serum diluted with PBS, contains a low concentration of IgG antibodies against *Borrelia burgdorferi*. Addition of 0.01 % methylisothiazolone and 0.01 % bromonitrodioxane. Ready-to-use.

**5.4. Calibrator C (Weak Positive Control)**

2 mL, human serum diluted with PBS, contains a medium concentration of IgG antibodies against *Borrelia burgdorferi*. Addition of 0.01 % methylisothiazolone and 0.01 % bromonitrodioxane. Ready-to-use.

**5.5. Calibrator D (Positive Control)**

2 mL, human serum diluted with PBS, contains a high concentration of IgG antibodies against *Borrelia burgdorferi*. Addition of 0.01 % methylisothiazolone and 0.01 % bromonitrodioxane. Ready-to-use.

**5.6. Enzyme Conjugate**

15 mL, anti-human-IgG-HRP (rabbit), in protein-containing buffer solution. Addition of 0.01 % methylisothiazolone and 0.01 % bromonitrodioxane and 5 mg/L Proclin™. Ready-to-use.

**5.7. Substrate**

15 mL, TMB (tetramethylbenzidine). Ready-to-use.

**5.8. Stop Solution**

15 mL, 1 N acidic solution. Ready-to-use.

**5.9. Sample Diluent**

60 mL, PBS/BSA buffer. Addition of 0.095 % sodium azide. Ready-to-use.

**5.10. Washing Buffer**

60 mL, PBS + Tween 20, 10x concentrate. Final concentration: dilute 1+9 with distilled water. If during the cold storage crystals precipitate, the concentrate should be warmed up at 37°C for 15 minutes.

**6. MATERIALS REQUIRED BUT NOT PROVIDED**

- 5 µL-, 100 µL- and 500 µL micro- and multichannel pipets
- Microtiter Plate Reader (450 nm)
- Microtiter Plate Washer
- Reagent tubes for serum dilution
- Deionized water
- Plastic Bag

## 7. SAMPLE COLLECTION AND HANDLING

Principally serum or plasma (EDTA, heparin) can be used for the determination. Serum is separated from the blood, which is aseptically drawn by venipuncture, after clotting and centrifugation. The serum or plasma samples can be stored refrigerated (2-8°C) for up to 7 days. For a longer storage they should be kept at -20°C. The samples should not be frozen and thawed repeatedly. Lipemic, hemolytic or bacterially contaminated samples can cause false positive or false negative results.

For the performance of the test the samples (not the standards) have to be diluted 1:101 with ready-to-use sample diluent (e.g. 5 µL serum + 500 µL sample diluent).

## 8. ASSAY PROCEDURE

### 8.1. Preparation of Reagents

**Washing Solution:** dilute before use 1+9 with deionized water. If during the cold storage crystals precipitate, the concentrate should be warmed up at 37°C for 15 minutes.

- Strict adherence to the protocol is advised for reliable performance. Any changes or modifications are the responsibility of the user.
- All reagents and samples must be brought to room temperature before use, but should not be left at this temperature longer than necessary.
- A standard curve should be established with each assay.
- Return the unused microtiter strips to the plastic bag and store them dry at 2-8°C.

### 8.2. Assay Steps

1. Prepare a sufficient amount of microtiter wells for the standards, controls and samples as well as for a substrate blank.
2. Pipet 100 µL each of the **diluted** (1:101) samples and the **ready-to-use** standards and controls respectively into the wells. Leave one well empty for the substrate blank.
3. Cover plate with the re-usable plate cover and incubate at room temperature for 60 minutes.
4. Empty the wells of the plate (dump or aspirate) and add 300 µL of diluted washing solution. This procedure is repeated totally three times. Residuals of the washing buffer are afterwards removed by gentle tapping of the microtiter plate on a tissue cloth.
5. Pipet 100 µL each of ready-to-use conjugate into the wells. Leave one well empty for the substrate blank.
6. Cover plate with the re-usable plate cover and incubate at room temperature for 30 minutes.
7. Empty the wells of the plate (dump or aspirate) and add 300 µL of diluted washing solution. This procedure is repeated totally three times. Residuals of the washing buffer are afterwards removed by gentle tapping of the microtiter plate on a tissue cloth.
8. Pipet 100 µL each of the ready-to-use substrate into the wells. This time also the substrate blank is pipetted.
9. Cover plate with the re-usable plate cover and incubate at room temperature for 20 minutes in the dark (e.g. drawer).
10. To terminate the substrate reaction, pipet 100 µL each of the ready-to-use stop solution into the wells. Pipet also the substrate blank.
11. After thorough mixing and wiping the bottom of the plate, perform the reading of the absorption at 450 nm (optionally reference wavelength of 620 nm). The color is stable for at least 60 minutes.

## 9. EVALUATION

### 9.1. Evaluation (Cut-off)

The calculated absorptions for the sera, as mentioned above, are compared with the value for the cut-off standard. If the value of the sample is higher, there is a positive result. For a value below the cut-off standard, there is a negative result. It seems reasonable to define a range of  $\pm 20\%$  around the value of the cut-off as a grey zone. In such a case the repetition of the test with the same serum or with a new sample of the same subject, taken after 2-4 weeks, is recommended. Both samples should be measured in parallel in the same run.

The positive control must show at least the double absorption compared with the cut-off standard.

### 9.2. Evaluation (U/mL)

The ready-to-use standards and controls of the Borrelia IgG antibody kit are defined and expressed in arbitrary units (U/mL). Consequently for a given subject follow-up controls become possible. The values for controls and standards in units are printed on the QC data sheet.

For this evaluation the absorptions of the standards and controls are graphically drawn *point-to-point* against their concentrations. From the resulting reference curve the concentration values for each sample can then be extracted in relation to their absorptions. It is also possible to use automatic computer programs. As curve fit *point-to-point* has to be chosen.

Calibrator B with its concentration of 10 U/mL serves as cut-off standard. Analogous to the cut-off evaluation a range of  $\pm 20\%$  around the cut-off is defined as a grey zone. Thus results between 8 and 12 U/mL are reported as borderline.

## 10. ASSAY CHARACTERISTICS

Borrelia burgdorferi ELISA	IgG
Intra-Assay-Precision	4.8 %
Inter-Assay-Precision	7.7 %
Inter-Lot-Precision	2.8 – 10.6 %
Analytical Sensitivity	0.99 U/mL
Recovery	79 – 88 %
Linearity	77 – 126 %
Cross-Reactivity	No cross-reactivity to TBEV (FSME)
Interferences	No interferences to bilirubin up to 0.3 mg/mL, hemoglobin up to 8.0 mg/mL and triglycerides up to 5.0 mg/mL
Clinical Specificity	97 %
Clinical Sensitivity	88 %

## 11. REFERENCES




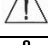

1. Ma B, Christen B, Leung D, Vigo-Pelfrey C: Serodiagnosis of Lyme borreliosis by immunoblot: Reactivity of various significant antibodies against *Borrelia burgdorferi*, J Clin Microbiol, 1992, 30:370-376
2. Münchhoff P, Wilske B, Preac-Mursic V, Schierz G: Antibodies against *Borrelia burgdorferi* in Bavarian forest workers, Zbl Bakt Hyg A, 1986, 263: 412-419
3. Dressler F, Whalen JA, Reinhardt BN, Steere AC: Western Blotting in the Serodiagnosis of Lyme Disease, J Infect Dis, 1993; 167:392-400
4. Tilton RC, Ryan RW: The Laboratory Diagnosis of Lyme Disease, J of Clin Immunoassay, 1993; 16:208-214
5. Barbour AG, Hayes SF: Biology of *Borrelia* species, Microbiol Rev, 1986, 50:381-400
6. Baranton G, Postic D, Saint Girons I, Boerlin P, Piffaretti JC, Assous M, Grimont PAD: Delineation of *Borrelia burgdorferi* sensu stricto, *Borrelia garinii* sp.nov. and group VS461 associated with Lyme borreliosis, Int J Syst Bacteriol 1992, 42:378-383
7. Berardi VP, Weeks KE, Steere AC: Serodiagnosis of early Lyme disease: Analyses of IgM and IgM antibody responses by using an antibody capture enzyme immunoassay, J Infect Dis, 1988, 158:754-760
8. Wilske B., Fingerle V., Herzer P., Hofmann A., Lehnert G., Peters H., Pfister H.-W., Preac-Mursic V., Soutschek E., Weber K.: Recombinant immunoblot in the serodiagnosis of Lyme Borreliosis, Med Microbiol Immunol, 1994, 182:255-270
9. Wilske B, Preac-Mursic V, Schierz G, Busch KV: Immunochemical and immunological analysis of European *Borrelia burgdorferi* strains, Zbl Bakt Hyg A, 1986, 263:92-102
10. Steere AC: Lyme Disease, New Engl J Med; 1989, 321:586-597
11. Grodzicki RL, Steere AC: Comparison of Immunoblotting and indirect enzyme-linked immunosorbent assay using different antigen preparations for diagnosing early Lyme disease, J Infect Dis, 1988, 157:790-797
12. Magnarelli LA, Anderson JF, Johnson RC: Cross-reactivity in serological tests for Lyme Disease and other spirochetal infections, J Infect Dis, 1987, 156:183-186

### Manufactured for:

Immuno-Biological Laboratories, Inc. (IBL-America)  
8201 Central Ave. NE, Suite P, Minneapolis, Minnesota 55432, USA  
Phone: +1 (763) - 780-2955 Fax.: +1 (763) - 780-2988  
Email: info@ibl-america.com Web: www.ibl-america.com



## SYMBOLS USED WITH IBL-AMERICA ASSAYS

Symbol	English	Deutsch	Francais	Espanol	Italiano
	European Conformity	CE-Konformitäts-kennzeichnung	Conforme aux normes européennes	Conformidad europea	Conformità europea
	Consult instructions for use	Gebrauchsanweisung beachten	Consulter les instructions d'utilisation	Consulte las Instrucciones	Consultare le istruzioni per l'uso
	In vitro diagnostic device	In-vitro-Diagnostikum	Usage Diagnostic in vitro	Diagnóstico in vitro	Per uso Diagnostica in vitro
	For research use only	Nur für Forschungszwecke	Seulement dans le cadre de recherches	Sólo para uso en investigación	Solo a scopo di ricerca
	Catalogue number	Katalog-Nr.	Référence	Número de catálogo	No. di Cat.
	Lot. No. / Batch code	Chargen-Nr.	No. de lot	Número de lote	Lotto no
	Contains sufficient for <n> tests/	Ausreichend für "n" Ansätze	Contenu suffisant pour "n" tests	Contenido suficiente para <n> ensayos	Contenuto sufficiente per "n" saggi
	Note warnings and precautions	Warnhinweise und Vorsichtsmaßnahmen beachten	Avertissements et mesures de précaution font attention	Tiene en cuenta y advertencias precauciones	Annoti avvisi e le precauzioni
	Storage Temperature	Lagerungstemperatur	Temperature conservation de	Temperatura conservacion de	Temperatura di conservazione
	Expiration Date	Mindesthaltbarkeits-datum	Date limite d'utilisation	Fecha de caducidad	Data di scadenza
	Legal Manufacturer	Hersteller	Fabricant	Fabricante	Fabbricante
<i>Distributed by</i>	Distributor	Vertreiber	Distributeur	Distribuidor	Distributore