



User's Manual

Mouse Interleukin-13(IL-13) ELISA

Sandwich Enzyme-Linked Immunosorbent Assay for the detection
of Mouse IL-13

REF

IB99558



96

Storage: 2-8°C

RUO

For research use only, not for use in diagnostic procedures.

Contents

1. Introduction	3
2. Principles of Method	3
3. Intended Use	4
4. Storage and Stability	4
5. Chemical Hazard	4
6. Kit Contents	5
7. Materials Required But Not Provided	6
8. Reagent preparation	6
1) Mouse IL-13 standard	6
2) Secondary Antibody	6
3) Streptavidin- HRP	7
4) Washing buffer	7
9. Assay Procedure	8
10. Characteristics	9
1) Typical result	9
2) Sensitivity	10
3) Specificity	11
4) Precision	11
5) Recovery	11
11. Troubleshooting	12
12. Reference	13

1. Introduction

Interleukin 13 (IL-13) is a protein that in humans is encoded by the IL13 gene. IL-13 was first cloned in 1993 and is located on chromosome 5q31 with a length of 1.4kb. IL-13 and IL-4 exhibit a 30% of sequence similarity and have a similar structure. IL-13 is cytokine secreted by many cell types, but especially T helper type 2 (Th2) cells, that is a mediator of allergic inflammation and disease. IL-13 has effects on immune cells that are similar to those of the closely related cytokine IL-4. However, IL-13 is suspected to be a more central mediator of the physiologic changes induced by allergic inflammation in many tissues. Although IL-13 is associated primarily with the induction of airway disease, it also has anti-inflammatory properties. IL-13 induces a class of protein-degrading enzymes, known as matrix metalloproteinases (MMPs), in the airways. These enzymes are required to induce egression of effete parenchymal inflammatory cells into the airway lumen where they are then cleared. Among other factors, IL-13 induces these MMPs as part of a mechanism that protects against excessive allergic inflammation that predisposes to asphyxiation. IL-13 is known to induce changes in hematopoietic cells, but these effects are probably less important than that of IL-4. Furthermore, IL-13 can induce immunoglobulin E (IgE) secretion from activated human B cells. Interestingly, deletion of IL-13 from mice does not markedly affect either Th2 cell development or antigen-specific IgE responses induced by potent allergens. In comparison, deletion of IL-4 abrogates these responses. Thus, rather than a lymphoid cytokine, IL-13 acts more prominently as a molecular bridge linking allergic inflammatory cells to the non-immune cells in contact with them, thereby altering physiological function. IL-13 induces its effects through a multi-subunit receptor that includes the alpha chain of the IL-4 receptor (IL-4R α) and at least one of two known IL-13-specific binding chains.

2. Principles of Method

The design of this assay is based on a sandwich Enzyme-Linked Immunosorbent Assay (ELISA). The microtiter plate provided in this kit has been pre-coated with a monoclonal antibody specific to Mouse IL-13. Samples are pipetted into these wells. Nonbound Mouse IL-13 and other components of the sample should be removed by washing, then monoclonal antibody specific to Mouse IL-13 added. In order to quantitatively determine the amount of Mouse IL-13 present in the sample, Streptavidin Horseradish Peroxidase (HRP) should be added to each microplate well. The final step, a TMB-substrate solution added to each well. Finally, a sulfuric acid solution is added and the resulting yellow colored product is measured at 450nm. Since the increases in absorbency is directly proportional to the amount of captured Mouse IL-13.

3. Intended Use

The IBL-AMERICA Mouse IL-13 ELISA kit is to be used for the determination of Mouse IL-13 in mice serum, mice plasma, cell lysate, culture supernatants and buffered solution. The assay will recognize native and recombinant Mouse IL-13.

This kit has been configured for research use only and is not to be used in diagnostic procedures.

4. Storage and Stability

All kit components of this kit are stable at 2 to 8°C. Any unused reconstituted standard should be discarded or frozen at -70°C. Standard can be frozen and thawed one time only without loss of immunoreactivity.

5. Chemical Hazard

- Stop solution: This reagent is an irritant to eyes, skin and mucous membranes. Avoid contact with eyes, skin and clothing. Wear suitable protective clothing, gloves and eye protection. In the event of contact with eyes or skin, wash immediately with plenty of water.
- Standard protein and 2nd Antibody containing Sodium Azide as a preservative.

6. Kit Contents

Contents	Number	Volume
96 Well Plate	1 (in aluminum foil bag with desiccant)	
Washing Buffer	2	(20X) 25 ml
Standard Protein	1 Glass vial (lyophilized)	
Standard/Sample Dilution Buffer	1	25 ml
Secondary Antibody	1 Glass vial (lyophilized)	
Streptavidin HRP(X100)	1	150 µl
Secondary antibody/ Streptavidin HRP Dilution Buffer	1	25 ml
Substrate (TMB)	1	15 ml
Stop Solution	1	15 ml
Protocol booklet	1	
Plate sealers	2	

① 96 Well Plate

: Mouse IL-13 microtiter plate, one plate of 96 wells (8 well strips x 12).

A plate using break-apart strips coated with a mouse monoclonal antibody specific to Mouse IL-13.

② Standard Protein

: Recombinant Mouse IL-13.

③ Secondary Antibody

: Biotinylated anti Mouse IL-13.

④ Streptavidin HRP

: Streptavidin Horseradish Peroxidase (HRP, enzyme)

⑤ Substrate (Stabilized chromogen)

: Tetramethylbenzidine (TMB) solution

⑥ Stop Solution

: 1N solution of sulfuric acid (H₂SO₄).

⑦ Plate sealer

: Adhesive sheet.

- Do not mix or interchange different reagents from various kit lots.

7. Materials Required But Not Provided

- ① Microtiter plate reader capable of measurement at or near 450nm.
- ② Calibrated, adjustable precision pipettes, preferably with disposable plastic tips (A manifold multi-channel pipette is desirable for large assays.)
- ③ Distilled or deionized water
- ④ Data analysis and graphing software
- ⑤ Vortex mixer
- ⑥ Polypropylene tubes for diluting and aliquoting standard
- ⑦ Absorbent paper towels
- ⑧ Calibrated beakers and graduated cylinders of various sizes

8. Reagent Preparation

1) Mouse IL-13 standard

1. Reconstitute the lyophilized Mouse IL-13 standard by adding 1 ml of *Standard/Sample Dilution Buffer* to make the 2 ng/ml standard stock solution. Allow solution to sit at RT for 5 minutes, then gently vortex to mix completely. Use within one hour of reconstituting
2. Prepare 1 ml of 500 pg/ml top standard by adding 250 µl of the above stock solution in 750 µl of *Standard/Sample Dilution Buffer*. Perform 2-fold serial dilutions of the top standards to make the standard curve within the range of this assay (7.8125 pg/ml ~ 500 pg/ml) as below. *Standard/Sample Dilution Buffer* serves as the zero standard (0 pg/ml).

Standard	Add	Into
500 pg/ml	250.0 µl of the std.(2 ng/ml)	750.0 µl of the Standard/Sample Dilution Buffer
250 pg/ml	500 µl of the std.(500 pg/ml)	500.0 µl of the Standard/Sample Dilution Buffer
125 pg/ml	500 µl of the std.(250 pg/ml)	500.0 µl of the Standard/Sample Dilution Buffer
62.5 pg/ml	500 µl of the std.(125 pg/ml)	500.0 µl of the Standard/Sample Dilution Buffer
31.25 pg/ml	500 µl of the std.(62.5 pg/ml)	500.0 µl of the Standard/Sample Dilution Buffer
15.625 pg/ml	500 µl of the std(31.25 pg/ml)	500.0 µl of the Standard/Sample Dilution Buffer
7.8125 pg/ml	500µl of the std.(15.625 pg/ml)	500.0 µl of the Standard/Sample Dilution Buffer
0 pg/ml	1.0 ml of the Standard/Sample Dilution Buffer	

2) Secondary Antibody

100X secondary antibody solution can be made by adding 150 µl *Secondary antibody/Streptavidin HRP dilution buffer* in the vial.

1. Equilibrate to room temperature, mix gently.

2. Mix 20 μ l *Secondary Antibody concentrated solution* (100X) + 1.98 ml *Secondary antibody/Streptavidin HRP dilution buffer*. (Sufficient for two 8-well strip, prepare more if necessary)

Label as “Working Secondary antibody Solution”.

3. Return the unused *Secondary Antibody concentrated solution* to the refrigerator.

3) Streptavidin HRP(X100)

1. Equilibrate to room temperature, mix gently.
2. Mix 20 μ l *Streptavidin HRP concentrated solution* (100X) + 1.98 ml *Secondary antibody/Streptavidin HRP dilution buffer*. (Sufficient for two 8-well strip, prepare more if needed)

Label as “Working Streptavidin HRP Solution”.

3. Return the unused *Streptavidin HRP concentrated solution* to the refrigerator.

4) Washing buffer

1. Equilibrate to room temperature, mix to re-dissolve any precipitated salt.
2. Mix 0.5 volume *Wash buffer concentrate solution* (20X) + 9.5 volumes of deionized water. Label as “Working Washing Solution”.
3. Store both the concentrated and the Working Washing Solution in the refrigerator.

* Directions for washing

1. Fill the wells with 300 μ l of “Working Washing Buffer”.
Let soak for 1 to 3 minutes and then all residual wash-liquid must be drained from the wells by aspiration (taking care not to scratch the inside of the well) or decantation, followed by forceful tapping of the plate on absorbent paper. Never insert absorbent paper directly into the wells.
If using an automated washer, the operating instructions for washing equipment should be carefully followed.
2. Incomplete washing will adversely affects the assay and renders false results.
3. It is recommended to use laboratory tape to hold plate strips to the plate frame while performing the plate washing to avoid strips coming free of the frame.

5) Sample preparation

Blood should be collected by venipuncture. For plasma samples, blood may be drawn into tubes containing sodium citrate or heparin, EDTA. The serum or plasma should be separated from the coagulated or packed cells by centrifugation. Samples may be shipped at room temperature and then stored refrigerated at 2-8°C if testing is to take place within one week after collection. If testing is to take place later than one week, samples should be stored at -20°C. Avoid repeated freeze/thawing.

9. Assay Procedure

- Allow all reagents to reach room temperature before use. Gently mix all liquid reagents prior to use.
 - All standards, controls and samples should be run in duplicate for confirmation of reproducibility.
 - A standard curve must be run with each assay.
 - If particulate matter is present in the analyte, centrifuge or filter prior to analysis.
 - Maintain a consistent order of components and reagents addition from well to well. This ensures equal incubation times for all wells.
- 1) Determine the number of 16-well strips needed for assay. Insert these in the flame(s) for current use (Re-bag extra strips and frame. Refrigerate for further use).
 - 2) For the standard curve, add 100 μ l of the standard to the appropriate microtiter wells. Add 100 μ l of the *Standard/Sample/secondary antibody Dilution Buffer* to zero wells.
 - 3) Serum and plasma require at least 20 fold dilution in the *Standard/Sample Dilution Buffer*. And add 100 μ l of samples to each well.
 - 4) Cover the plate with the plate cover and incubate for 2 hours at 37°C.
 - 5) Thoroughly aspirate or decant the solution from the wells. Wash the wells 3 times (See “Directions for washing”).
 - 6) Pipette 100 μ l of “Working Secondary Antibody Solution” into each well.
 - 7) Cover the plate with the plate cover and incubate for 1 hour at 37°C.
 - 8) Thoroughly aspirate or decant the solution from the wells. Wash the wells 3 times (See “Directions for washing”).
 - 9) Add 100 μ l “Working Streptavidin HRP Solution” to each well.
 - 10) Cover the plate with the plate cover and incubate for 30 minutes at 37°C.
 - 11) Thoroughly aspirate or decant the solution from the wells. Wash the wells 3 times (See “Directions for washing”).
 - 12) Add 100 μ l of *Substrate* to each well. The liquid in the wells should begin to turn blue.
 - 13) Incubate the plate at room temperature.
 - Do not cover the plate with aluminum foil, or color may develop.
The incubation time for chromogen substrate is often determined by the microtiter plate reader used. O.D. values should be monitored and the substrate reaction stopped before O.D. of the positive wells exceeds the limits of the instrument. O.D. values at 450nm can only be read after the Stop Solution has been added to each well.
 - Because the *Substrate* is light sensitive, avoid the remained *Substrate* solution prolonged exposure to light.

- Typically, reaction is stopped 5~10 minutes after treatment of Substrate, but this time can be adjusted as the user desires.
- 14) Add 100 µl of *Stop Solution* to each well. The solution in the wells should change from blue to yellow.
 - 15) Read the absorbance of each well at 450nm. Read the plate within 20 minutes of adding the *Stop Solution*.
 - 16) Plot on graph paper the absorbance of the standard against the standard concentration (Optimally, the background absorbance can be subtracted from all data points, including standards, unknowns and controls, prior to plotting.). Draw a smooth curve through these points to construct the standard curve.
 - 17) Read the Mouse IL-13 concentrations for the unknown samples and controls from the standard curve plotted in step 18. Multiply value(s) obtained for the unknown sample by the dilution factor (Samples producing signals greater than that of the highest standard should be further diluted in the *Standard/Sample Dilution Buffer*).

10. Characteristics

1) Typical result


The standard curve below is for illustration only and **should not be used** to calculate results in your assay.

A standard curve must be run with each assay.

Standard Mouse IL-13 (pg/ml)	Optical Density (at 450nm)
0	0.109
7.8125	0.154
15.615	0.193
32.15	0.267
62.5	0.418
125	0.684
250	1.121
500	1.742

< Limitations >

- Do not extrapolate the standard curve beyond the 500 pg/ml standard point.
- Other buffers and matrices have not been investigated.
- The rate of degradation of native Mouse IL-13 in various matrices has not been investigated.

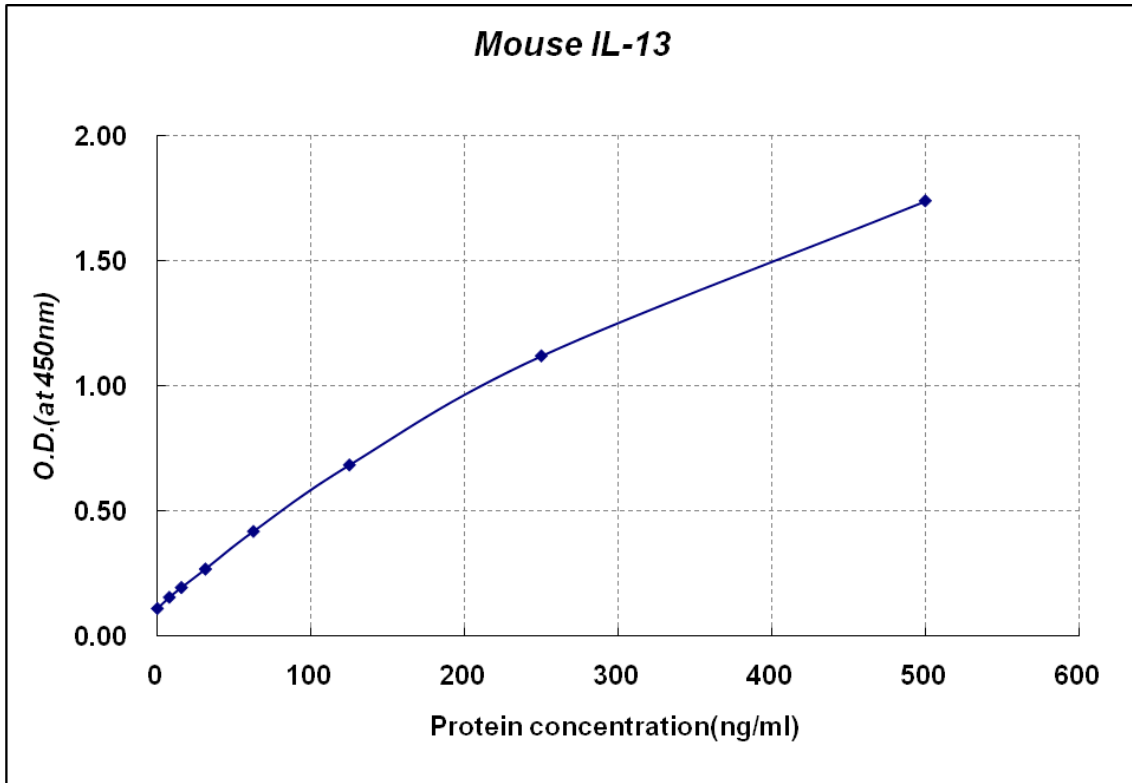

(TMB reaction incubate at room temperature for 5 min)

2) Sensitivity

The minimal detectable dose of Mouse IL-13 was calculated to be 7.153 pg/ml, by subtracting three standard deviations from the mean of 12 zero standard replicates (ELISA buffer, S0) and intersecting this value with the standard curve obtained in the same calculation.

N	1	2	3	4	5	6	7	8	9	10	11	12
ZERO	0.115	0.097	0.093	0.089	0.095	0.089	0.089	0.090	0.091	0.101	0.098	0.094

Average	SD	LLD	LLD mean(pg/ml)
0.095	0.007	0.117	7.153

3) Specificity

The following substances were tested and found to have no cross-reactivity: Human IL-13, Rat IL-13.

4) Precision

① Within-Run (Intra-Assay)

(n=12)

Mean (pg/ml)	SD (pg/ml)	CV (%)
64.298	4.484	6.974
123.888	6.148	4.963
250.395	9.333	3.727
500.439	32.929	6.580

② Between-Run (Inter-Assay)

(n=12)

Mean (pg/ml)	SD (pg/ml)	CV (%)
58.220	2.643	4.540
129.262	11.782	9.115
248.314	21.587	8.693
501.537	38.727	7.722

5) Recovery

Recovery on addition is 91.668~95.390% (mean 93.077%)

Added Analyte (pg/ml)	Serum(1/20)+added analyte (450nm)	Serum(1/20)(450nm) +added analyte(450nm)	Recovery (%)
61.527	53.265	49.114	92.2072
128.377	103.692	95.052	91.668
247.771	225.535	209.842	93.042
500.648	447.864	427.218	95.390

11. Troubleshooting

Problem	Possible Cause	Solution
High signal and background in all wells	<ul style="list-style-type: none"> • Insufficient washing 	<ul style="list-style-type: none"> • Increase number of washes • Increase time of soaking between in wash
	<ul style="list-style-type: none"> • Too much Streptavidin-HRP 	<ul style="list-style-type: none"> • Check dilution, titration
	<ul style="list-style-type: none"> • Incubation time too long 	<ul style="list-style-type: none"> • Reduce incubation time
	<ul style="list-style-type: none"> • Development time too long 	<ul style="list-style-type: none"> • Decrease the incubation time before the stop solution is added
No signal	<ul style="list-style-type: none"> • Reagent added in incorrect order, or incorrectly prepared 	<ul style="list-style-type: none"> • Review protocol
	<ul style="list-style-type: none"> • Standard has gone bad (If there is a signal in the sample wells) 	<ul style="list-style-type: none"> • Check the condition of stored standard
	<ul style="list-style-type: none"> • Assay was conducted from an incorrect starting point 	<ul style="list-style-type: none"> • Reagents allows to come to 20~30°C before performing assay
Too much signal – whole plate turned uniformly blue	<ul style="list-style-type: none"> • Insufficient washing –unbound AV-HRP remaining 	<ul style="list-style-type: none"> • Increase number of washes carefully
	<ul style="list-style-type: none"> • Too much Streptavidin -HRP 	<ul style="list-style-type: none"> • Check dilution
	<ul style="list-style-type: none"> • Plate sealer or reservoir reused, resulting in presence of residual Streptavidin -HRP 	<ul style="list-style-type: none"> • Use fresh plate sealer and reagent reservoir for each step
Standard curve achieved but poor discrimination between point	<ul style="list-style-type: none"> • Plate not developed long enough 	<ul style="list-style-type: none"> • Increase substrate solution incubation time
	<ul style="list-style-type: none"> • Improper calculation of standard curve dilution 	<ul style="list-style-type: none"> • Check dilution, make new standard curve
No signal when a signal is expected, but standard curve looks fine	<ul style="list-style-type: none"> • Sample matrix is masking detection 	<ul style="list-style-type: none"> • More diluted sample recommended
Samples are reading too high, but standard curve is fine	<ul style="list-style-type: none"> • Samples contain protein levels above assay range 	<ul style="list-style-type: none"> • Dilute samples and run again
Edge effect	<ul style="list-style-type: none"> • Uneven temperature around work surface 	<ul style="list-style-type: none"> • Avoid incubating plate in areas where environmental conditions vary • Use plate sealer

12. Reference

- 1) Minty A, Chalon P, Derocq JM, Dumont X, Guillemot JC, Kaghad M, Labit C, Leplatois P, Liauzun P, Miloux B (March 1993). "Interleukin-13 is a new human lymphokine regulating inflammatory and immune responses". *Nature* 362 (6417): 248–50.
- 2) McKenzie AN, Culpepper JA, de Waal Malefyt R, Brière F, Punnonen J, Aversa G, Sato A, Dang W, Cocks BG, Menon S (April 1993). "Interleukin 13, a T-cell-derived cytokine that regulates human monocyte and B-cell function". *Proc. Natl. Acad. Sci. U.S.A.* 90 (8): 3735–9.
- 3) Morgan JG, Dolganov GM, Robbins SE, Hinton LM, Lovett M (October 1992). "The selective isolation of novel cDNAs encoded by the regions surrounding the human interleukin 4 and 5 genes". *Nucleic Acids Res.* 20 (19): 5173–9.
- 4) Zurawski G, de Vries JE (January 1994). "Interleukin 13, an interleukin 4-like cytokine that acts on monocytes and B cells, but not on T cells". *Immunol. Today* 15 (1): 19–26.
- 5) Wynn TA (2003). "IL-13 effector functions". *Annu. Rev. Immunol.* 21: 425–56.

◆ Ordering Information

For orders, please contact :

Immuno-Biological Laboratories, Inc. (IBL-America)

Address: 8201 Central Ave NE, Suite P, Minneapolis, MN 55432

Toll Free: (888) 523-1246

Fax: (763) 780-2988

E-mail: info@ibl-america.com

Website: www.ibl-america.com

For technical advice, please contact:

E-mail : info@ibl-america.com

Website : www.ibl-america.com