

Code No. 29020

**Anti-Human
MDA5 Rabbit IgG Affinity Purify**Volume : 100 μ g

Introduction : A full-length cDNA encoding MDA5 was isolated by screening a melanoma differentiation-subtracted library and a placenta cDNA library, also known as Ifih1(Interferon induced with helicase C domain protein 1). Yoneyama et al. have identified a DExD/H box-containing RNA helicase, retinoic acid inducible gene I (RIG-I), as an essential component of the sensor of intracellular dsRNA (1). Furthermore, two other DExD/H-box-containing RNA helicases, MDA5 and LGP2, were found in the mammalian database (2). MDA5 acts as a positive regulator in the virus-induced activation of type I IFN genes. On the other hand, LGP2 functions as a dominant-negative regulator of RIG-I/MDA5-mediated signaling. These molecules have critical roles for antiviral innate immune responses.

Antigen : Synthetic peptide of part of N terminal of human MDA5 (TIEDRNRIAAENNGN)

Purification : Purified with antigen peptide

Form : Lyophilized product from 1% BSA in PBS containing 0.05% NaN₃

How to use : 1.0 mL deionized water will be added to the product (the conc. comes up 100 μ g/mL)

Stability : Lyophilized product, 5 years at 2 – 8 °C
: Solution, 2 years at –20 °C


Application : Western blotting in concentration of 1 – 5 μ g/mL

Specificity : Reacts with human and mouse MDA5

Reference : 1. Yoneyama M, Kikuchi M, Natsukawa T, Shinobu N, Imaizumi T, Miyagishi M, Taira K, Akira S, Fujita T. The RNA helicase RIG-I has an essential function in double-stranded RNA-induced innate antiviral responses. *Nat Immunol.* 2004 Jul; 5 (7):730-7.
2. Yoneyama M, Kikuchi M, Matsumoto K, Imaizumi T, Miyagishi M, Taira K, Foy E, Loo YM, Gale M Jr, Akira S, Yonehara S, Kato A, Fujita T. Shared and unique functions of the DExD/H-box helicases RIG-I, MDA5, and LGP2 in antiviral innate immunity. *J Immunol.* 2005 Sep 1;175 (5):2851-8.

Example : lane 1:
HeLa cell lysate without IFN β stimulation
lane 2:
HeLa cell lysate with IFN β stimulation

When this antibody is used in Western Blotting, several bands appear although the antibody specifically detects the MDA5 protein. We recommend to perform positive and negative controls in parallel.

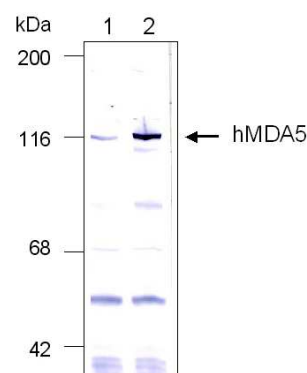

For research use only, not for use in diagnostic procedures.

Distributed by:



Immuno-Biological Laboratories, Inc. Toll-Free: 888-523-1246
8201 Central Ave NE, Suite P Email: info@IBL-America.com
Minneapolis, MN 55432 Web: www.IBL-America.com